

# Detect

## Detect Covid-19 Test™ Covid-19 Molecular Home Test



## Instructions for Use For Healthcare Providers

For in vitro diagnostic use.

For use under Emergency Use Authorization only.

# Contents

<b>Intended Use</b>	<b>2</b>
<b>Summary and Explanation of the Test</b>	<b>3</b>
<b>Principles of the Procedure</b>	<b>4</b>
<b>Assay/Reagents</b>	<b>5</b>
<b>Warnings and Precautions</b>	<b>6</b>
<b>Operating Conditions</b>	<b>9</b>
<b>Procedure</b>	<b>10</b>
<b>Quality Control</b>	<b>31</b>
<b>Retests</b>	<b>31</b>
<b>Limitations</b>	<b>32</b>
<b>Performance Characteristics</b>	<b>34</b>
<b>Bibliography</b>	<b>42</b>
<b>Symbols and Abbreviations</b>	<b>42</b>

# 1. Intended Use

The Detect Covid-19 Test™ (the Detect test) is a molecular *in vitro* diagnostic test for the qualitative detection of nucleic acid from the novel coronavirus SARS-CoV-2 that causes Covid-19.

This test is authorized for non-prescription home use with self-collected anterior nasal (nasal) swab samples from individuals aged 14 years or older suspected of Covid-19. This test is also authorized for non-prescription home use with adult-collected anterior nasal swab samples from individuals aged 2 years or older suspected of Covid-19.

This test is also authorized for non-prescription home use with self-collected anterior nasal (nasal) swab samples from individuals aged 14 years or older, or adult collected anterior nasal swab samples from individuals aged 2 years or older, without symptoms or other epidemiological reasons to suspect Covid-19 when tested twice over three days with at least 24 hours (and no more than 48 hours) between tests.

SARS-CoV-2 viral RNA is generally detectable in anterior nasal swab specimens during the acute phase of infection. Positive results indicate the presence of viral RNA, but clinical correlation with past medical history and other diagnostic information is necessary to determine infection status. Positive results do not rule out bacterial infection or co-infection with other viruses. Individuals who test positive with the Detect Covid-19 Test should self-isolate and seek follow-up care with their physician or healthcare provider as additional testing may be necessary.

Negative results should be treated as presumptive and may be confirmed with a molecular assay performed in a laboratory, if necessary for patient management. Negative results do not preclude SARS-CoV-2 infection and should not be used as the sole basis for treatment or management decisions for the individual, including infection control decisions. Negative results should be considered in the context of an individual's recent exposures, history and the presence of clinical signs and symptoms consistent with Covid-19.

Individuals who test negative and continue to experience Covid-19 like symptoms of fever, cough and/or shortness of breath may still have SARS-CoV-2 infection and should seek follow up care with their physician or healthcare provider.

Test results will be reported to relevant public health authorities in accordance with local, state, and federal requirements, using appropriate LOINC and SNOMED codes, as defined by the Laboratory In Vitro Diagnostics (LIVD) Test Code Mapping for SARS-CoV-2 Tests provided by the CDC. Automatic test result reporting will be performed by the Detect App™ and the Detect secure cloud server.

The Detect Covid-19 Test is authorized for non-prescription self-test by individuals aged 14 years or older and/or, as applicable, for an adult lay user testing another person aged 2 years or older in a non-laboratory setting. The Detect Covid-19 Test is only for use under the Food and Drug Administration's Emergency Use Authorization.

## **2. Summary and Explanation of the Test**

An outbreak of pneumonia of unknown etiology in Wuhan City, Hubei Province, China was initially reported to the World Health Organization (WHO) in December 2019. Chinese authorities identified a novel coronavirus SARS-CoV-2 (cause of Covid-19 respiratory disease) which has resulted in confirmed human infections worldwide, including the United States. Cases of severe respiratory illness and deaths have been reported. Patients can become infected with SARS-CoV-2 virus through contact with a contaminated environment or person.

The Detect Covid-19 Test is a molecular in vitro diagnostic test that aids in the diagnosis of Covid-19 through the identification of the SARS-CoV-2 RNA in nasal swab specimens.

In asymptomatic individuals (those without Covid-19 symptoms), the Detect Covid-19 Test should be used as a serial test.

### **What is serial testing?**

Serial testing involves testing the same person multiple times within a few days. Such testing for Covid-19 increases the chances of identifying infections earlier and should be used for people who are not exhibiting any symptoms.

### **How do I use the Detect Covid-19 Test for serial testing?**

If you're asymptomatic and not exposed to Covid-19 and your first test is negative, you take a second test after at least 24 hours but within 48 hours. If your first test is positive, then you are likely to have Covid-19 currently and should consult with a healthcare provider without waiting to use the second test. If only your second test is positive, then you are also likely to have Covid-19 and should consult with a healthcare provider.

### 3. Principles of the Procedure

The Detect test uses RT-LAMP (Reverse Transcription Loop-mediated Isothermal Amplification) and lateral flow strip technologies to recognize nucleic acids from the Open Reading Frame 1ab (ORF1ab) region of the SARS-CoV-2 genome. The Detect test also identifies nucleic acids from a human gene that serves as a control for sample collection, extraction, reagent integrity, and test execution.

Isothermal amplification occurs at elevated temperature within a disposable tube placed into the reusable Detect Hub™. After amplification, the tube is inserted into the Reader and the tube's liquid wicks onto the lateral flow strip. On the lateral flow strip's sample pad, SARS-CoV-2 and control amplicons bind colored particles and flow through the lateral flow strip's membrane, where they are captured by immobilized antibodies at distinct lines on the strip. A valid negative result must show the Sample Processing Control line. A positive result must show the SARS-CoV-2 line and may or may not also show the Sample Processing Control line.

In asymptomatic patients, serial testing is required to assist in identifying infected individuals and facilitate timely infection control practices. A negative test result does not rule out infection but repeat testing done on a weekly basis may decrease the risks of false negative results.

For asymptomatic users, an initial negative test result should be the first of a minimum of two tests. An asymptomatic individual undergoing serial testing with two or more negative results may require ongoing serial testing or confirmatory testing with this or a different SARS-CoV-2 test, depending on patient history and potential exposures. An asymptomatic individual undergoing serial testing with one or more positive results indicates that SARS-CoV-2 RNA is present but does not rule out coinfection with other pathogens.


## 4. Assay/Reagents

### 4.1 Materials

The Detect Covid-19 Test contains enough reagents to process one self-collected sample.


#### Materials Provided

- Swab (sterilized)
- Test Tube (contains Collection Buffer)
- Detect Cap™ (contains lyophilized reagent bead)
- Dropper (contains buffer)
- Reader (contains lateral flow strip inside of plastic housing)

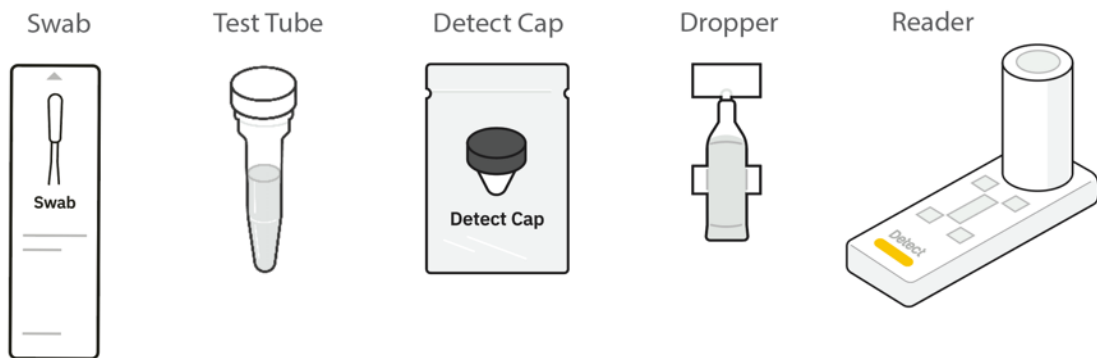

#### Required but Not Provided

- Detect Hub (Model 21101, device sold separately)
- Test instructions from one of the following:
  - Detect App (free) and a smartphone - visit [detect.com/app](https://detect.com/app) to see a list of compatible devices and download the app
  - Detect Covid-19 Test Instructions for Use - visit [detect.com/ifu](https://detect.com/ifu) to download

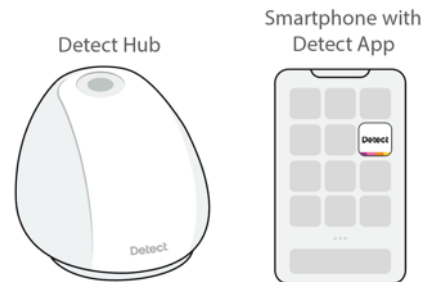

#### Minimum smartphone requirements for the Detect App:

iOS Models released after 2018 and using the operating system iOS 13 or higher.

Android Models released after 2018 using the Android API level 29 (version 10 or higher) including models from Samsung, Google, and OnePlus. The Detect App is not compatible with Windows smartphones, Android tablets, or iPads.

## 5. Warnings and Precautions

### 5.1 General



- For *in vitro* diagnostic use.
  - For use under FDA Emergency Use Authorization only.
  - This product has not been FDA cleared or approved but has been authorized by FDA under an EUA.
  - This product has been authorized only for the detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens.
  - The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of Covid-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated, or authorization is revoked sooner.
  - For more information on EUAs please visit:  
<https://www.fda.gov/emergencypreparednessand-response/mcm-legal-regulatory-and-policy-framework/emergency-useauthorization>
  - For the most up to date information on Covid-19, please visit:  
<https://www.cdc.gov/coronavirus/2019-ncov/index.html>
  - The intensity of the test line does not necessarily correlate to the amount of SARS-CoV-2 virus in the sample.
  - Samples should be tested as quickly as possible after sample collection.
  - Failure to follow the instructions for use may adversely affect test performance and/or invalidate the test result.
  - Only use the test components provided. Do not use swabs from other tests. Do not re-use any of the components included in the Detect Covid-19 Test. Only the Detect Hub may be re-used for testing other samples.
  - Positive results are indicative of the presence of SARS-CoV-2 RNA.
- Treat all biological specimens, including used test components, as if capable of transmitting infectious agents. Because it is often impossible to know which might be infectious, all biological specimens should be handled using standard precautions. Guidelines for specimen handling

are available from the U.S. Centers for Disease Control and Prevention [1, 2] and the Clinical and Laboratory Standards Institute [3].

- Do not ingest.
- Keep out of reach of children.
- Avoid contact with skin and eyes.
- Do not apply the Test Tube buffer directly onto the skin or mucous membranes or ingest. If contact with the body occurs, rinse with water. If irritation persists, seek medical advice. If swallowed: call a poison center/doctor if you feel unwell.
- Do not apply the buffer contained in the Dropper directly onto the skin or mucous membranes or ingest. If contact with the body occurs, rinse with water. If irritation persists, seek medical advice. If swallowed: call a poison center/doctor if you feel unwell.

## 5.2 Storage & Handling

- Store all components at 59 °F to 86 °F (15 °C to 30 °C).
- Do not open components until you are ready to perform testing.
- Open all packages carefully to avoid losing small components.
- Do not use Test Tubes that are wet or have leaked or spilled.
- Do not use the Detect Cap if its storage pouch is punctured or not fully sealed. The cap contains a freeze-dried bead of reagents that is sensitive to moisture.
- Do not use the Detect Covid-19 Test past the Use By date on the test box label.
- All components other than the Detect Hub are single-use and should be disposed of after use.
- Follow all instructions carefully as shown in the Detect App or the Detect Covid-19 Test Instructions for Use.
- Touch only the plastic handle of the Swab with your hands to avoid contaminating the soft tip of the Swab.
- Do not insert the Swab deeper than 1-2 cm into your nose. A deeper swabbing will not yield more accurate results.



- See Section 8 for detailed instructions on the control included with the test. It helps indicate whether the reaction is taking place correctly.

### 5.3 Components & Reagents

- Do not remove the Detect Cap after screwing it onto the Test Tube.
- Do not use a Detect Cap that has been dropped after removing it from the packaging.
- Do not shake the Test Tube except as described in Section 7, Step 4-5 (Test Tube Preparation).
- Begin Test Tube processing for each sample within 1 hour of collecting the sample.
- Place swab immediately into the Test Tube after collecting samples. Failure to do so may result in dried swabs and yield an incorrect test result.
- Do not put anything into the chimney of the Reader until instructed to do so. Doing so may lead to indeterminate results.
- ② • Each single-use Swab is used for one test. Do not reuse Swabs or use a Swab other than the one provided in the test.
- ② • Each single-use Test Tube is used for one test. Do not reuse Test Tubes.
- ② • Each single-use Detect Cap is used for one test. Do not reuse Detect Caps.
- ② • Each single-use Dropper is used for one test. Do not reuse Droppers.
- ② • Each single-use Reader is used for one test. Do not reuse Readers.
- Do not process the samples using any protocol other than the one described in Section 7, Step 4-6 as other protocols have not been tested.
- Do not tamper with the Reader or attempt to remove the Test Tube once inserted into the Reader.

## 6. Operating Conditions

- The test should be used between 59 °F and 86 °F (15 °C and 30 °C). Failure to do so may yield invalid or inaccurate results.
- The test is best used in a room with adequate lighting and away from glare. Failure to do so may result in an inability to see the results on the test.
- The Hub must be run on a level surface and should not be moved during operation. Failure to do so may yield invalid or inaccurate results.
- If a power failure occurs or if the Hub is unplugged while the Test Tube is in the Hub, the test result is invalid and the user should be retested following the retest procedure described in Section 9.2.

## 7. Procedure

### Step 1: Obtain items required but not provided in the test

You will need the items listed below to run your Detect Covid-19 Test. These items are not included in the Detect test.

- Detect Hub. You can purchase the Hub from Detect, Inc. at [detect.com](https://detect.com).
- Smartphone. Go to [detect.com/app](https://detect.com/app) for the list of compatible smartphones.
- The Detect App installed on your smartphone. Download the Detect App from [detect.com/app](https://detect.com/app).

Instead of using a smartphone with the Detect App, you can also download the Detect Covid-19 Test Instructions for Use from [detect.com/ifu](https://detect.com/ifu). If this is your first time taking the Detect Covid-19 Test, using the Detect App is recommended.

### Step 2: Prepare to run your test

Carefully read the Detect Covid-19 Test Getting Started instructions (on the inner lid of the test box) before you run your test. These directions will help you download the Detect App and complete the test correctly and safely.

If you do not understand the instructions, do not run the test. Contact Detect customer support at [support@detect.com](mailto:support@detect.com) or call toll-free at 855-322-3692 for help.

### Step 3: Open the Detect App on your smartphone and follow the instructions

The Detect App will use pictures, videos, and on-screen instructions to guide you through collecting your nasal sample and running the Detect Covid-19 Test. Be sure to carefully follow all the instructions given within the Detect App. If you do not, your test may give an invalid result.

If you do not understand the directions, contact Detect customer support at [support@detect.com](mailto:support@detect.com) or call toll-free at 855-322-3692 for help.

- You may need to update your Detect App to the latest version before running your test. Follow any on-screen prompts to update the Detect App.
- Within the Detect App, tap on the “Start new test” button.

## **Step 4: Run the Detect test**

Follow on-screen instructions to run the Detect test.

### **Step 4-1: Get ready to run the Detect test**

Review and complete the instructions on the “Getting ready” screens:

- Check that the test has not expired by reviewing the Use By date on the box.
- Wipe down your work surface and wash your hands, being careful to rinse them thoroughly.
- Plug in the Detect Hub.
- Set aside time (approximately 65 minutes) to run the test.
- Tap “Next step” to continue to the test activation step.


### **Step 4-2: Activate the Detect test**

Open your test and find the activation code on the inner lid of the box. Use your smartphone camera to scan the activation code or manually enter the 7-word activation code printed on the inner lid of your test box.


### **Step 4-3: Gather the components required to run the test**

Place the Prepare Pack pouch and the divider on a flat surface. See the Detect App video showing materials that you need to run a test as shown [here](#).

- Unpack the Swab, Detect Cap, and Test Tube from the Prepare Pack pouch.

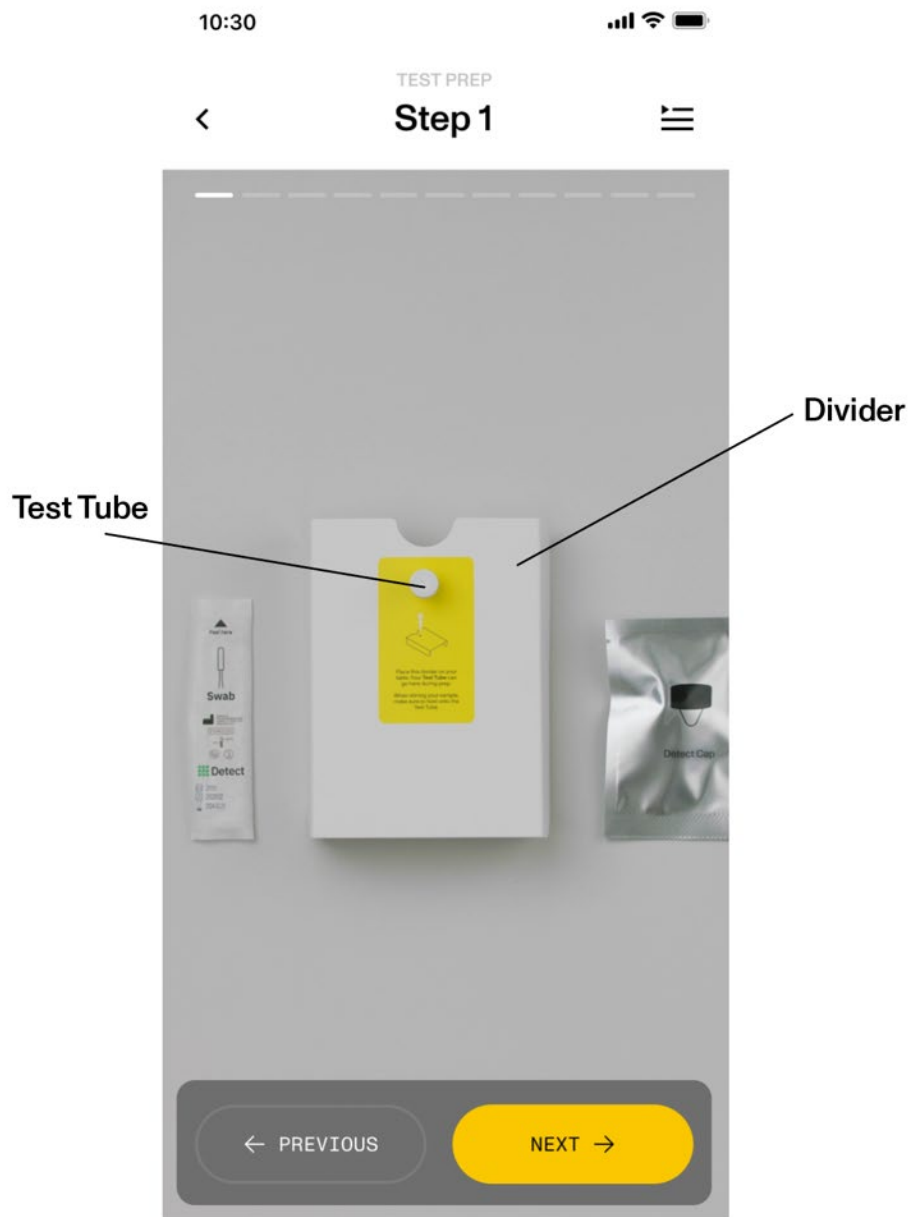

Follow the directions printed on the divider to place the Test Tube in the divider.

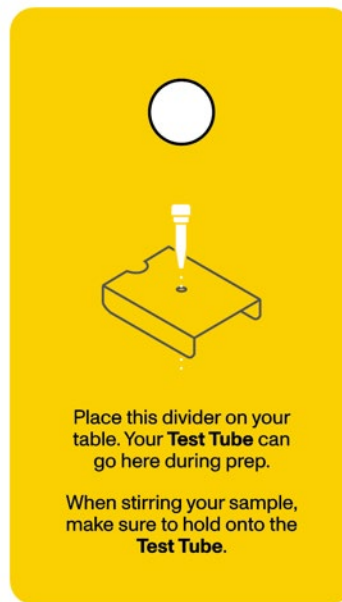

**REMINDER:** Do not open the Detect Cap pouch yet.

#### Step 4-4: Prepare and collect a nasal sample

- Unscrew the cap from the Test Tube. Be careful since the Test Tube contains liquid. Place the Test Tube back into the divider as shown in the Detect App video.




- Remove the Swab on the side that says “peel here”. Make sure that the soft tip of the Swab does not touch anything besides your nose.
- The Detect App video will demonstrate how to collect a nasal sample.

**REMINDER:** To properly collect a nasal sample, swab both of your nostrils with the same Swab following the directions below.


The Detect test is for use with self-collected samples for individuals aged 14 years and older and for samples collected and tested by an adult caregiver for individuals aged 2 years and older.

- For adult collection, insert the Swab into the nostril until just the soft tip is completely inside (about 1 inch). For pediatric collection insert the Swab into the child's nostril only ½ inch. Swab in a circle around the inside wall of the nostril 5 times. Then gently remove the Swab, insert into your other nostril, and repeat. Make sure the Swab stays in full contact with the inside of the nostril.

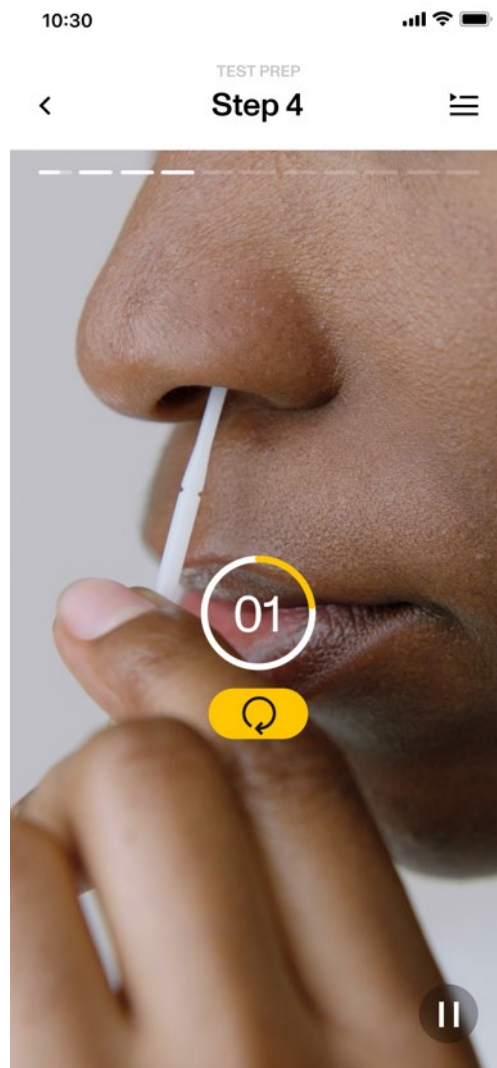



- Hold the Test Tube with one hand and fully submerge the Swab tip in the Test Tube and vigorously twirl the Swab for 15 seconds, as shown in the Detect App video. Do not break off the Swab tip. Discard the Swab by placing it back inside its package.

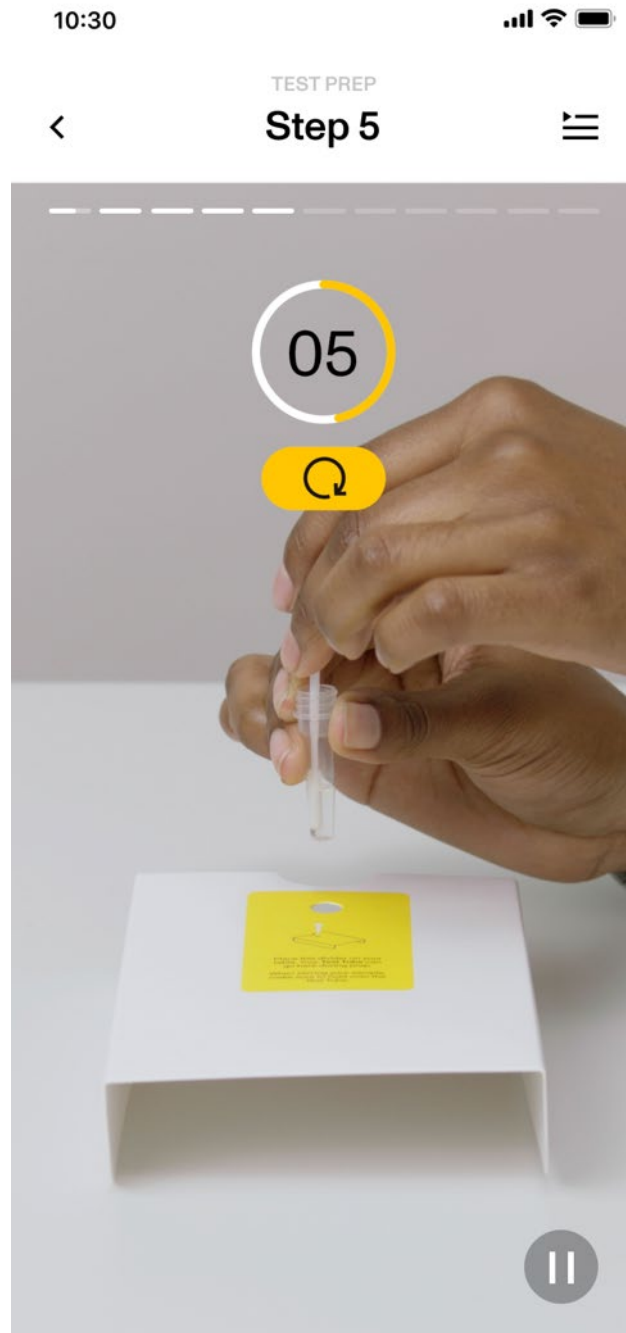


#### **Step 4-5: Test Tube Preparation**


**REMINDER:** Pick up the Detect Cap pouch. While it's closed, gently push the contents inside down toward the bottom of the pouch. The Detect Cap can fly out of the pouch if opened too quickly.

- Slowly and carefully open the pouch and gently remove the Detect Cap. You may need to remove the silica packet first.
- On the inside of the Detect Cap you'll see a small white reagent bead—this contains the reagents to run your test.
- Screw the Detect Cap onto the Test Tube, tightening as much as possible. Do not open the Test Tube after this step.
- Hold the Test Tube by the cap, turn it upside down, and shake it vigorously side to side for 10 seconds, as shown in the Detect App video. Then turn the tube right side up again, continuing to hold it by the cap.

**REMINDER:** Do not shake the Test Tube up and down. Shake it side to side to ensure the reagent bead contained within the Detect Cap is always in contact with the liquid. This will help dissolve the reagents needed to run the Detect test.

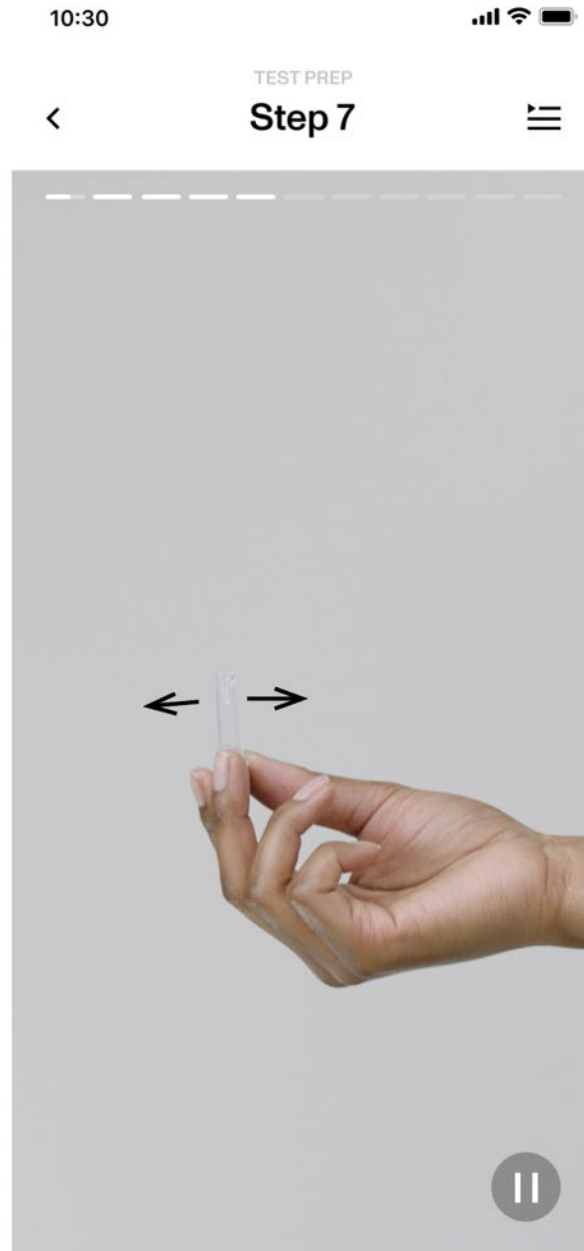


- Forcefully bring your whole arm downward to move the liquid to the bottom of the tube. Do this a few times to ensure the liquid is not stuck in the cap or against the walls of the tube.

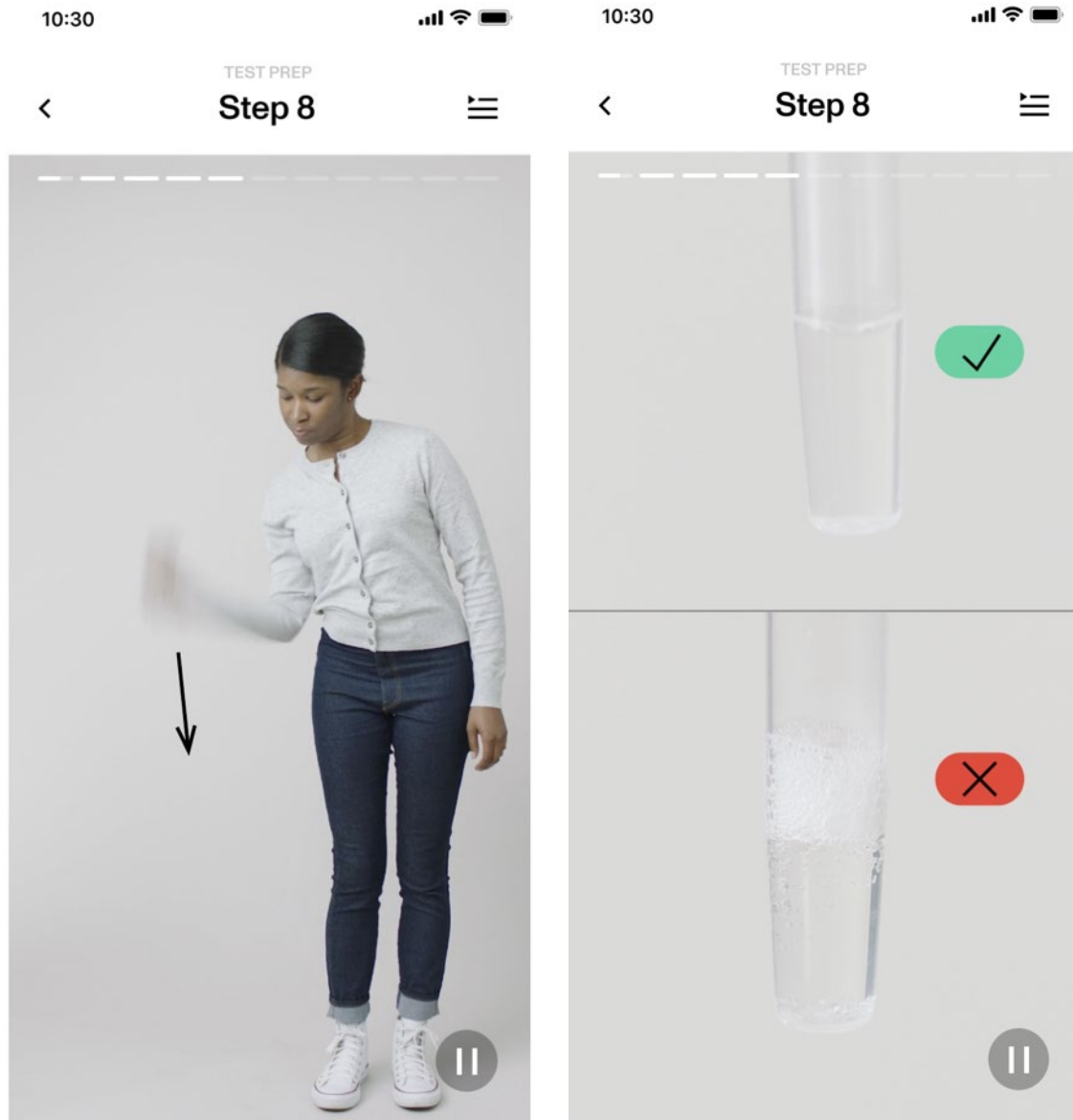

**REMINDER:** Liquid clinging to tube walls or cap may result in an invalid test.

#### Step 4-6: Process the sample

- Place the Test Tube into the well of the Detect Hub, pushing all the way in, as shown in the Detect App video. The Hub will beep once and the green light will start blinking. Your sample will start processing automatically and will take 55 minutes. The blinking green light will turn solid when it's complete.

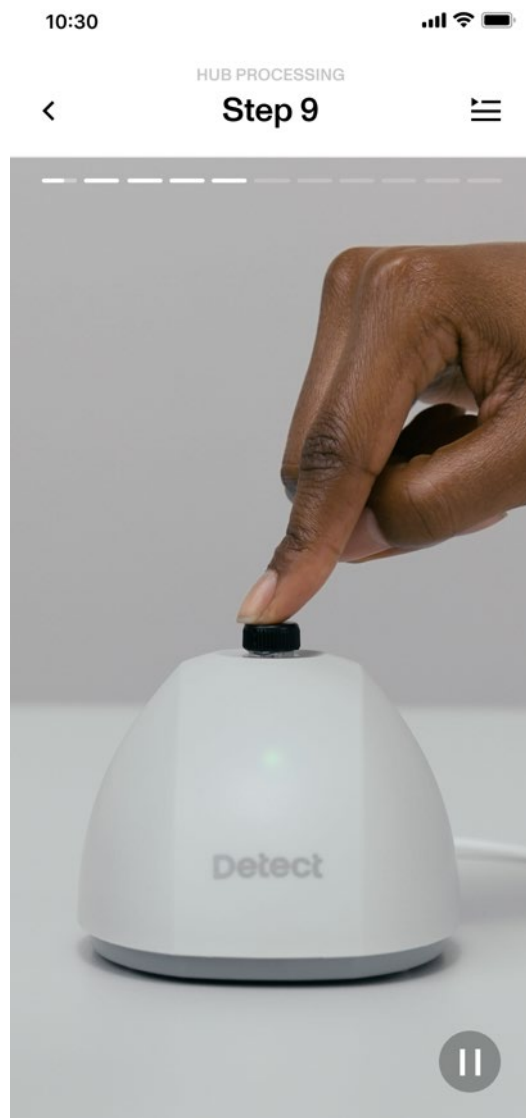

**REMINDER:** Do not remove the Test Tube from the Hub until it is time to use it.

#### Step 4-7: Results

- Pick up the Dropper by the rectangular tip and snap your wrist downward to collect the liquid at the bottom as shown in the Detect App video. Carefully twist off the tip of the Dropper. Hold the Dropper gently to avoid squeezing it and spilling the liquid.
- Insert the dropper as far as it will go into the Reader chimney to avoid spilling, and using both hands, squeeze it firmly to dispense all the liquid, as shown in the Detect App video. Discard the Dropper.



- Remove the Test Tube from the Detect Hub, and place it in the chimney of the reader. Using both thumbs, press down on the tube, as shown in the Detect App video. You may hear a pop as the blade at the bottom of the chimney opens the tube, allowing the liquid inside to flow through the Reader.

10:30




TEST RESULTS



Step 14



**REMINDER:** Keep pressing until the tube is completely flush with the lip of the Reader chimney. Liquid should begin to flow through the Reader and should be visible within the Reader's window in about 10 seconds. If there's no flow, firmly tap the Reader against a hard surface 3 times to help trigger flow.

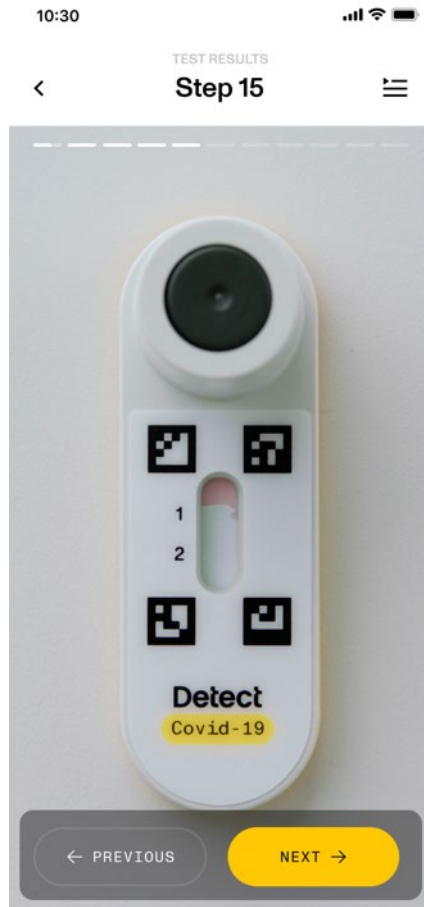


- The red lines will take approximately 10 minutes to fully develop. Line darkness and fullness may vary.

**REMINDER: Do not tamper with the Reader or try to remove the Test Tube from the Reader.**

- The Detect App will then guide you through identifying the visible lines on the Reader to interpret the result of the test.
- Line 1 is the test line, which tells you whether SARS-CoV-2 (the virus that caused Covid-19) was identified in your sample.
- Line 2 is the control line, which tells you if the test was performed correctly for a negative test result.



- As shown below on the right, lines can be faint or incomplete.

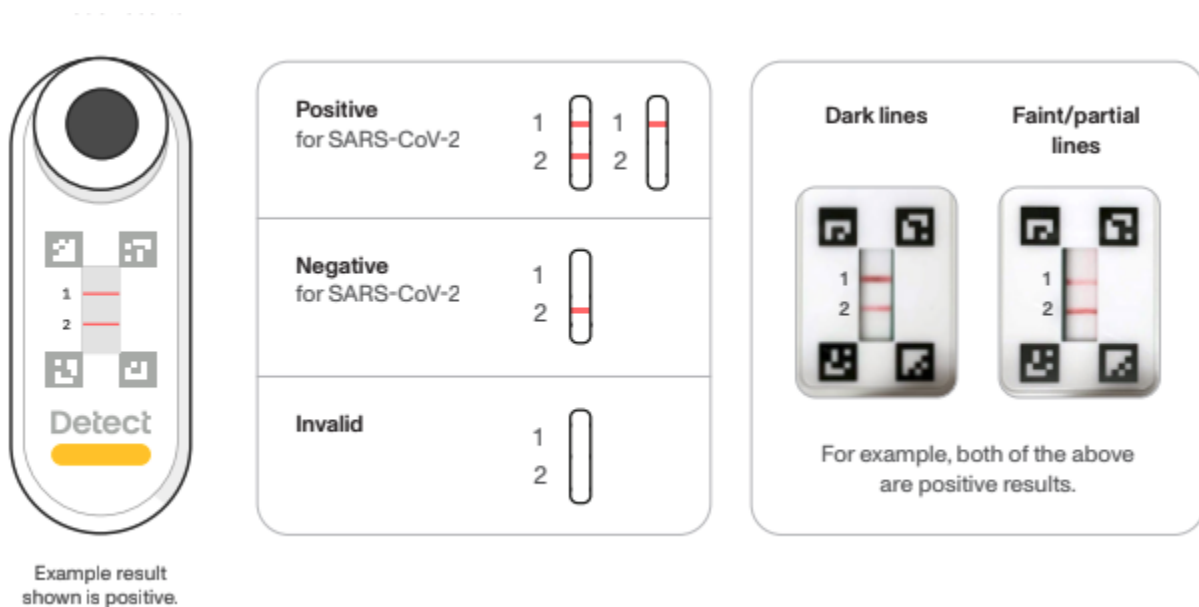


### Step 5: Understanding the test result

The app will show results as Negative, Positive, or Invalid. Please see Detect's Fact Sheet for Individuals ([resources.detect.com](https://resources.detect.com)) for more information on understanding the test result.

- Regardless of the test result, it is important that when you are sick you practice social distancing and good hygiene.
- If you develop symptoms or your symptoms persist or become more severe, if you are concerned about your health, or if you develop one of the emergency warning signs ([www.cdc.gov/coronavirus](https://www.cdc.gov/coronavirus)), then you should seek medical attention immediately.

## Step 5-1: Understanding a Negative result

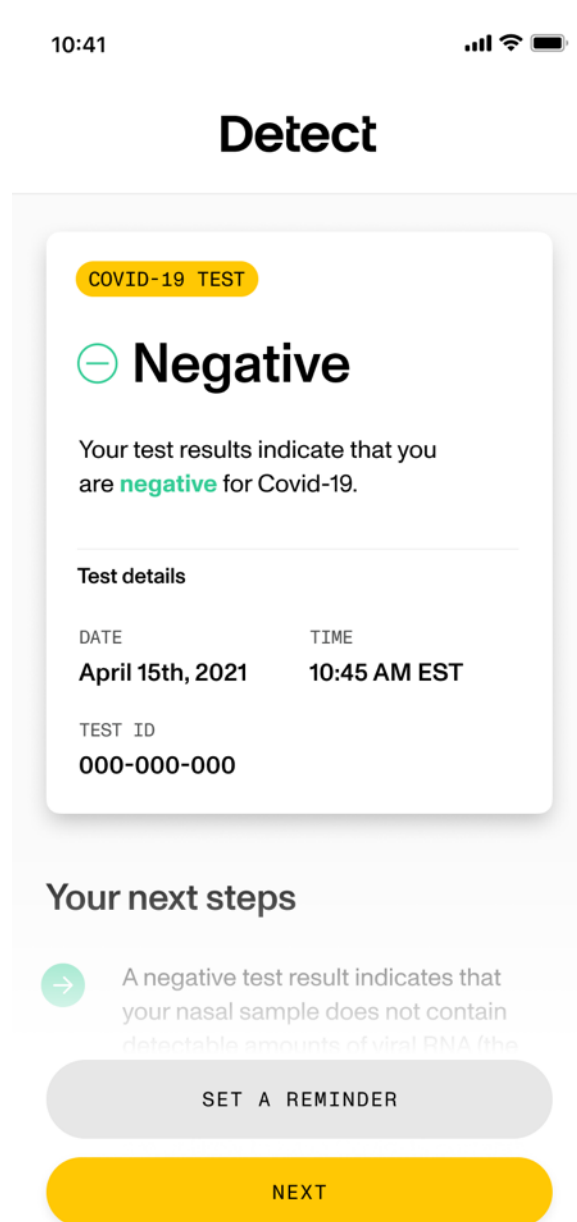




A Negative result means that the Detect Covid-19 Test did not identify the SARS-CoV-2 virus that causes Covid-19 in your sample and it is unlikely that you currently have a Covid-19 infection.

- If you are symptomatic and your symptoms persist or become more severe, and you obtained a negative result for your sample, seek help from your HCP even if the results remain negative based on serial testing.


- If you have no symptoms or reasons to suspect Covid-19 infection, this was your first result, and your result was negative, you must collect and test a second sample at least 24 hours after your first test and within the next 48 hours.
- In the following scenarios it is possible that the test may give a negative result even if you have Covid-19 (called a false negative result).
  - A false negative can occur if the sample was not collected or processed properly.
  - if you are too early or late in your Covid-19 infection to accurately identify a low amount of SARS-CoV-2 virus in your sample.
  - A false negative can also occur if the SARS-CoV-2 virus's genetic material changes (mutates) such that the Detect Covid-19 Test cannot identify the virus.
- For serial testing programs, additional testing is required when a negative result is obtained for the first sample. Additional testing may also be necessary if the individual was exposed to someone who tested positive for SARS-CoV-2 (the virus that can cause Covid-19), or in communities with high numbers of positive cases (high prevalence of infection).
- Even if you do not have Covid-19, you may still have another type of infection. Many other viruses can cause similar symptoms to Covid-19, and these may be the cause of your current illness.
- Tell your healthcare provider whether or not you have symptoms.

## Step 5-2: Understanding a Positive result

10:41   

# Detect

COVID-19 TEST



## Positive


Your test results indicate that you are **positive** for Covid-19.

### Test details

DATE	TIME
April 15th, 2021	10:45 AM EST

TEST ID  
000-000-000

### Your next steps

 A positive test result indicates that your nasal sample contains detectable amounts of viral RNA (the genomic material from the virus that causes Covid-19). This means you are


NEXT

A Positive result means the Detect Covid-19 Test identified the SARS-CoV-2 virus that causes Covid-19 in your sample, and it is very likely that you currently have a Covid-19 infection.

- You should self-isolate at home per CDC recommendations to stop spreading the virus to others. Consult the CDC recommendations regarding self-isolation at [cdc.gov/coronavirus](https://www.cdc.gov/coronavirus).
- Consult your healthcare provider as soon as possible and tell him or her that you tested positive for Covid-19 using the Detect Covid-19 Test.

- Tell your healthcare provider if you have symptoms or no symptoms.
- Tell your healthcare provider to view the Detect Covid-19 Test Fact Sheet for Healthcare Professionals at [resources.detect.com](https://resources.detect.com).
- There is a small possibility that this test can give a positive result that is wrong (a false positive result) particularly when used in a population without many cases of COVID-19 infection. Your healthcare provider will work with you to determine how best to care for you based on the test results along with medical history, and your symptoms
- There is still a chance of co-infection with another type of illness.
- If you do not have any symptoms, particularly if you live in an area with low numbers of Covid-19 infections and have had no exposure to anyone diagnosed with Covid-19, additional testing to confirm your result may be required.

## Step 5-3: Understanding an Invalid result

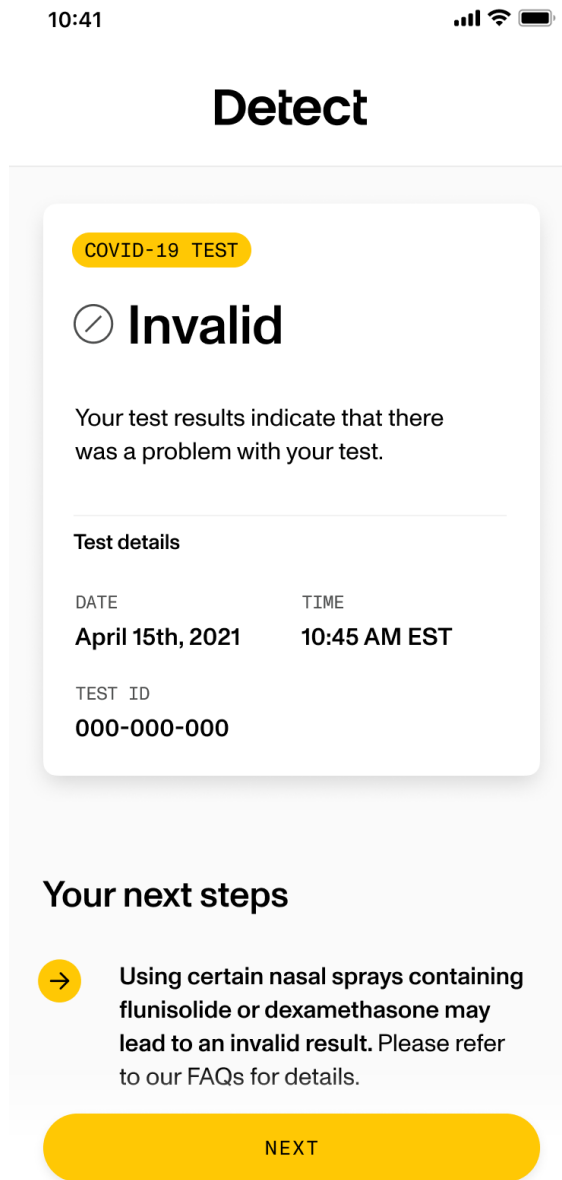

An Invalid result means that there was an error, and the Detect Covid-19 Test was unable to provide a result. You will need to perform a retest.

Some common causes of Invalid results are:

- The sample was improperly collected and did not include enough nasal material. Make sure to thoroughly swab both nostrils as directed in Step 4-4.

- The Test Tube was not prepared correctly (either the reagent was not properly dissolved or the liquid was not fully collected into the bottom of the Test Tube). See Step 4-5 for details.
- The processing step was run incorrectly. Make sure you fully insert the Test Tube into the Detect Hub until it beeps and do not remove the Test Tube until the processing step is complete. See Step 4-6 for details. A power failure or interruption during test processing can also cause an invalid result.
- If the result is invalid, retest. You must use a new Detect Covid-19 Test and a new Swab. Contact Detect customer support at [support@detect.com](mailto:support@detect.com) or call toll-free at 855-322-3692 for a replacement test.

## 8. Quality Control

### 8.1 Internal Control

#### **CONTROL**

Each reaction includes a Sample Processing Control (SPC) that is shown as a separate line in the reader.

**Sample Processing Control (SPC, Line 2)** - Ensures that the sample was processed correctly. The Sample Processing Control is designed to amplify a human control gene that will be present in a swab sample collected by a subject. The SPC determines whether sample collection was performed correctly and whether amplification reaction conditions were appropriate (temperature, time, and reagent mixing). The SPC should be positive in a human sample that tests negative and positive in a human sample that tests positive; however, a positive sample that lacks the SPC is still valid.

## 9. Retests

### 9.1 Reasons to Retest

If an **INVALID** test result occurs, repeat the test once according to instructions in 9.2 (Retest Procedure). If the repeat test fails to produce a valid result, please contact Detect Customer Support at 855-322-3692 or [support@detect.com](mailto:support@detect.com).

If you have no reason to suspect Covid-19 infection, this was your first result, and your result was negative, you must collect and test a second sample at least 24 hours after your first test and within the next 48 hours.

### 9.2 Retest Procedure

1. Prepare the work surface as detailed in Section 7.1 (Setup).
2. Obtain a new Detect Covid-19 Test.
3. Repeat the test procedure outlined in Section 7 (Procedure).



## 10. Limitations

- Performance has only been evaluated for self-collected nasal swab samples. Use of the Detect Covid-19 Test with other specimen types (such as saliva) has not been evaluated.
- A false negative result may occur if the individual's nose is swabbed incorrectly. False negative results may also occur if the amount of SARS-CoV-2 virus present on the swab is below the test's limit of detection.
- False negative results may occur in patients currently taking high dose biotin (vitamin B7) supplements. Biotin levels of 0.88 µg/mL or higher in the nasal sample may result in incorrect test results (false negatives).
- Invalid results may occur if:
  - The patient is currently using certain nasal sprays used to fight allergies, such as dexamethasone or flunisolide-containing products. Dexamethasone and flunisolide levels higher than 0.25 mg/mL and 7.5% volume, respectively, in the nasal sample may result in invalid test results in some cases.
  - High amounts of moisturizing lotions or hand soap comes in direct contact with the swab and is introduced into the test reaction. Hands should therefore be washed and well-rinsed before performing this test.
- The performance of this test was established based on the evaluation of a limited number of clinical specimens. Clinical performance has not been established with all circulating variants but is anticipated to be reflective of the common variants in circulation at the time and location of the clinical evaluation. Performance at the time of testing may vary depending on the variants circulating, including newly emerging strains of SARS-CoV-2 and their prevalence, which change over time.
- This test cannot rule out diseases caused by other bacterial or viral pathogens.
- Test results should be interpreted together with the patient's medical history, clinical signs and symptoms, and the results of other diagnostic tests performed.
- As with other tests, negative results do not rule out SARS-CoV-2 infections and should not be used as the sole basis for patient management decisions.
- This is a qualitative test. Test line intensity is not related to the quantity of virus in the sample.

- The material identified by this test—viral nucleic acid—may persist in the body, even after the patient is not infectious anymore. A positive result on this test does not imply that the patient is infectious, or that the virus is causing the subject's clinical symptoms.

## 11. Performance Characteristics

### 11.1 Analytical Sensitivity (Limit of Detection)

The Limit of Detection (LoD) is the lowest amount of virus that can be determined in 95% of samples and was determined by testing the Detect test's analytical sensitivity with heat-inactivated SARS-CoV-2 BEI, USA-WA1/2020) in pooled nasal matrix. Over two lots of tests, the LoD of the Detect Covid-19 Test was determined to be 313 copies per swab, which is equivalent to 800 copies per mL if all virus is transferred from the swab to the buffer.

Viral Load (genomic copies/swab)	Lot 1 SARS- CoV-2 Detection Rate	% Detected	Lot 2 SARS- CoV-2 Detection Rate	% Detected
313	20/20	100%	20/20	100%
156	20/20	100%	18/20	90%
78	15/20	75%	14/20	70%

### 11.2 Analytical Reactivity (Inclusivity)

An *in silico* inclusivity study was performed to analyze the Detect test's primer binding sequences in the SARS-CoV-2 genome to demonstrate that the primers will identify all variants of the SARS-CoV-2 virus identified to date (July 2021) and predict inclusivity of the Detect Covid-19 Test. A total of 2,397,589 sequences from the GISAID EpiCoV database ([www.gisaid.org](http://www.gisaid.org)) were evaluated in the study.

Based on the *in silico* analysis combined with laboratory testing, at least 98.6% of U.S. viral genomes from the past 90 days are expected to be robustly identified by the Detect test's SARS-CoV-2 primer set.

### 11.3 Analytical Specificity/Exclusivity (Cross-Reactivity)

The Detect test's cross-reactivity with closely related pathogens, common disease agents, and normal and pathogenic flora that may be present in the respiratory tract was tested in triplicate by spiking the organism or genomic material from the organism directly into Detect reactions at the concentrations listed in the table below. The Detect test showed no interaction with any of the 31 organisms tested. The Detect test was also tested repeatedly with pooled human nasal matrix (>500 replicates in all) and showed no cross-reactivity.

Organism	Target	Concentration Tested (in final reaction)	SARS-CoV-2 #Positive /# tested	Cross-reactivity with Detect
Human coronavirus 229E	Virus	1.00E+05 TCID <sub>50</sub> /mL	0/3	No
Human coronavirus OC43	Virus	1.00E+05 TCID <sub>50</sub> /mL	0/3	No
Human coronavirus HKU1	Synthetic RNA	6.85E+05 copies/mL	0/3	No
Human coronavirus NL63	Virus	4.00E+04 TCID <sub>50</sub> /mL	0/3	No
MERS-coronavirus	Virus	5.00E+03 TCID <sub>50</sub> /mL	0/3	No
SARS-coronavirus	Virus	1.50E+03 TCID <sub>50</sub> /mL	0/3	No
Adenovirus (Adenoid 71)	Virus	1.00E+05 TCID <sub>50</sub> /mL	0/3	No
Human Metapneumovirus (hMPV)	Virus	1.00E+05 TCID <sub>50</sub> /mL	0/3	No
Parainfluenza virus 1	Virus	5.00E+04 TCID <sub>50</sub> /mL	0/3	No
Parainfluenza virus 2	Virus	1.00E+05 TCID <sub>50</sub> /mL	0/3	No
Parainfluenza virus 3	Virus	1.00E+05 TCID <sub>50</sub> /mL	0/3	No

Parainfluenza virus 4	Virus	1.60E+04 TCID <sub>50</sub> /mL	0/3	No
Influenza A	Virus	1.00E+05 CEID <sub>50</sub> /mL	0/3	No
Influenza B	Virus	1.00E+05 TCID <sub>50</sub> /mL	0/3	No
Enterovirus 68	Virus	1.00E+05 TCID <sub>50</sub> /mL	0/3	No
Respiratory syncytial virus (Subgroup A)	Virus	1.00E+05 PFU/mL	0/3	No
Rhinovirus 89	Virus	8.00E+04 TCID <sub>50</sub> /mL	0/3	No
<i>Chlamydia pneumoniae</i>	Bacteria	1.00E+06 IFU/mL	0/3	No
<i>Haemophilus influenzae</i>	Bacteria	1.00E+06 CFU/mL	0/3	No
<i>Legionella pneumophila</i>	Bacteria	1.00E+06 CFU/mL	0/3	No
<i>Mycobacterium tuberculosis</i>	Bacteria	1.00E+06 CFU/mL	0/3	No
<i>Streptococcus pneumoniae</i>	Bacteria	8.80E+04 CFU/mL	0/9	No
<i>Streptococcus pyogenes</i>	Bacteria	1.00E+06 CFU/mL	0/3	No
<i>Bordetella pertussis</i>	Bacteria	1.00E+06 CFU/mL	0/3	No
<i>Mycoplasma pneumoniae</i>	Bacteria	1.00E+06 CFU/mL	0/3	No
<i>Pneumocystis jirovecii</i> (PJP)- <i>S. cerevisiae</i> *	Yeast	1.00E+06 CFU/mL	0/3	No
<i>Candida albicans</i>	Yeast	1.00E+06 CFU/mL	0/3	No
<i>Pseudomonas aeruginosa</i>	Bacteria	1.00E+06 CFU/mL	0/3	No
<i>Staphylococcus epidermis</i>	Bacteria	1.00E+06 CFU/mL	0/3	No

<i>Streptococcus salivarius</i>	Bacteria	1.00E+06 CFU/mL	0/3	No
---------------------------------	----------	--------------------	-----	----

\* Due to limited pathogen availability, cross-reactivity was tested with a recombinant version of *S. cerevisiae* containing genomic material from PJP.

In addition, *in silico* cross-reactivity analysis of Detect primer sequences was performed by comparing them to representative genomic sequences of the specific respiratory microorganisms below, downloaded from the NCBI database. The table below details all instances of ≥80% homology between a primer and respiratory microorganism genome.

Greater than 80% homology was only apparent for a single SARS-CoV-2 primer with *Pneumocystis jirovecii* (PJP) and two primers with *Candida albicans*. Further, none of the labelled primers required for identification of the amplified nucleic acid target showed ≥80% homology with any of the listed respiratory microorganism genomes. Therefore, *in silico* analysis identified no potential unintended cross-reactivity of the Detect test with the listed respiratory pathogens, including other coronaviruses.

Organism	SARS-CoV-2 primer set
Human coronavirus 229E	no alignment found
Human coronavirus OC43	no alignment found
Human coronavirus HKU1	no alignment found
Human coronavirus NL63	no alignment found
MERS-CoV	no alignment found
SARS-CoV	no alignment found
Adenovirus (e.g. C1 Ad. 71)	no alignment found
Human Metapneumovirus (hMPV)	no alignment found
Parainfluenza virus 1	no alignment found
Parainfluenza virus 2	no alignment found
Parainfluenza virus 3	no alignment found
Parainfluenza virus 4	no alignment found
Influenza A	no alignment found

Influenza B	no alignment found
Enterovirus (e.g. EV68)	no alignment found
Respiratory syncytial virus	no alignment found
Rhinovirus A	no alignment found
Rhinovirus B	no alignment found
Rhinovirus C	no alignment found
<i>Chlamydia pneumoniae</i>	no alignment found
<i>Haemophilus influenzae</i>	no alignment found
<i>Legionella pneumophila</i>	no alignment found
<i>Mycobacterium tuberculosis</i>	no alignment found
<i>Streptococcus pneumoniae</i>	no alignment found
<i>Streptococcus pyogenes</i>	no alignment found
<i>Bordetella pertussis</i>	no alignment found
<i>Mycoplasma pneumoniae</i>	no alignment found
<i>Pneumocystis jirovecii</i> (PJP)	single primer only, 83%
<i>Candida albicans</i>	Two primers at 82% and 89%
<i>Pseudomonas aeruginosa</i>	no alignment found
<i>Staphylococcus epidermidis</i>	no alignment found
<i>Staphylococcus salivarius</i>	no alignment found

## 11.4 Analytical Specificity (Interfering Substances)

Common endogenous and exogenous substances that might be present in clinical nasal swab samples were tested for interference with the Detect test. Each potentially interfering substance was spiked into both negative pooled nasal matrix and contrived positive pooled nasal matrix spiked with heat-inactivated SARS-CoV-2 virus at 2X LoD. From these pools, triplicate swabs were tested using the Detect test. The interfering substances and their concentrations are listed in the table below. The results show that the Detect test is robust to a wide range of potentially interfering substances.

Interfering Substance	Final Concentration in Nasal Matrix Pool	Negative Samples # Negative /# tested	Positive Samples # Positive /# tested	Interference Observed
Rhinocort Allergy	15% v/v	3/3	3/3	No
Afrin Nasal Congestion Relief Spray	15% v/v	3/3	3/3	No
Zicam Cold Remedy Nasal Spray	15% v/v	3/3	3/3	No
Chloraseptic Sore Throat Spray	15% v/v	3/3	3/3	No
Flonase Allergy Relief Nasal Spray	15% v/v	3/3	3/3	No
Mupirocin	1 mg/mL	3/3	3/3	No
Neo-Syneprine	15% v/v	3/3	3/3	No
Nasal Saline Spray	15% v/v	3/3	3/3	No
Tobramycin	600 µg/mL	3/3	3/3	No
Fresh whole blood	15%	3/3	3/3	No
Biotin	3.5 µg/mL	3/3	1/3*	Yes
	0.875 µg/mL	3/3	3/3	No
Dexamethasone	0.5 mg/mL	2/3**	3/3	Yes
	0.25 mg/mL	3/3	3/3	No
Flunisolide	15%	2/3**	3/3	Yes
	7.5% v/v	3/3	3/3	No
Mucin	1 mg/mL	3/3	3/3	No



Triamcinolone	15% v/v	3/3	3/3	No
Mometasone nasal spray	1 mg/mL	3/3	3/3	No
Method All-Purpose Surface Cleaner	15% v/v	3/3	3/3	No

\* The positive samples that could not be detected had false negative results.

\*\* The samples without a negative result had invalid results.

Other potentially interfering substances (hand soap and lotion) that might be present in the home environment were tested for interference under real-world use cases. Hand soap and lotion were tested by having users wash their hands with hand soap without rinsing or apply hand lotion immediately prior to grasping the swab by the tip. These intentionally contaminated swabs were then spiked with contrived positive nasal samples at 2x LoD and assayed on the Detect test. Neither of the tested substances interfered with the test's ability to identify SARS-CoV-2.

<b>Interfering Substance</b>	<b>Contamination delivery route</b>	<b>Positive Samples # Positive /# tested</b>	<b>Interference Observed</b>
Dial Antibacterial Liquid Hand Soap	On operator's hands	3/3	No
Aveeno Daily Moisturizing Lotion	On operator's hands	3/3	No

## 11.5 Clinical Evaluation

A prospective, multi-center clinical study was conducted in the United States in subject's homes or a simulated home environment. Testing was performed by the untrained subject or the subject's parent/guardian for children under 14. The study enrolled symptomatic subjects and asymptomatic subjects with a recent exposure, each of whom self-collected two nasal swab samples. For each subject, one swab was collected and sent to a reference laboratory and tested using a high sensitivity FDA-authorized SARS-CoV-2 RT-PCR test by trained laboratory personnel as a comparator, while the other swab was run through the Detect test by the untrained subject or parent/guardian. Comparing the Detect test's results to those produced by the high sensitivity FDA-authorized SARS-CoV-2 RT-PCR test, the Positive Percent Agreement (PPA) was 90.9.% (30/33) and the Negative Percent Agreement (NPA) was

97.5% (77/79). Both apparent false positives were due to subject's misinterpreting the test result; the Detect App has since been modified to significantly reduce the frequency of this user error.

### Clinical Study Results Summary

	FDA-authorized SARS-CoV-2 RT-PCR Assay	
Detect	Positive	Negative
Positive	30	2*
Negative	3**	77

\*Both of the apparent false positive results were due to subject misinterpretation using an older version of the Detect App. A modified version of the Detect App has since been developed and shown to be effective in reducing the frequency of subjects misinterpreting the test result.

\*\*One of the three apparent false negative samples gave a negative result when tested with a second highly sensitive EUA SARS-CoV-2 RT-PCR assay.

**Positive Percent Agreement (PPA):** 90.9% (30/33), (95% CI: 76.4%-96.9%)

**Negative Percent Agreement (NPA):** 97.5% (77/79), (95% CI: 91.2%-99.3%)

A second evaluation was conducted testing the ability of untrained lay users to correctly identify and interpret near-cutoff positive samples with the Detect test. Twelve subjects were each given three blinded contrived swab samples, a mix of negative and low-positive (1.9x LoD) to test. Comparing their results to those expected, the PPA was 94.4% (17/18) and the NPA was 100.0% (18/18), demonstrating that in the hands of lay users the Detect Covid-19 Test can reliably identify low-positive samples.

### Near-Cutoff Results Summary

	Contrived Sample Expected Result	
Detect	Positive	Negative
Positive	17	0
Negative	1	18

**Positive Percent Agreement (PPA):** 94.4% (17/18), (95% CI: 74.2%-99.0%)






**Negative Percent Agreement (NPA):** 100.0% (18/18), (95% CI: 82.4%-100.0%)

## 12. Bibliography


- 1 Centers for Disease Control and Prevention. "Interim Guidelines for Collecting, Handling, and Testing Clinical Specimens from Persons for Coronavirus Disease 2019 (Covid-19)." (Refer to latest edition.) <https://www.cdc.gov/coronavirus/2019-ncov/lab/guidelines-clinical-specimens.html>
- 2 Centers for Disease Control and Prevention. "Biosafety in Microbiological and Biomedical Laboratories (BMBL) 5th Edition." (Refer to latest edition.) <https://www.cdc.gov/labs/BMBL.html>
- 3 Clinical and Laboratory Standards Institute. "Protection of Laboratory Workers from Occupationally Acquired Infections; Approved Guideline." Document M29 <https://clsi.org/standards/products/microbiology/documents/m29/> (Refer to latest edition).

## 13. Symbols and Abbreviations

The following symbols are used throughout this manual:

Symbol	Definition
	Biohazard – Potentially infectious materials. Precautions must be observed.
	For <i>in vitro</i> diagnostic use
	Internal control
	Do not re-use
	Manufacturer

Manufacturer:

 Detect, Inc.  
530 Old Whitfield St.  
Guilford, CT 06437 USA

Customer Support: [support@detect.com](mailto:support@detect.com) or call toll-free +1 855-322-3692

---

**Detect**

# **Detect<sup>TM</sup> Hub User Manual**

**For in vitro diagnostic use.**

**For use with the Detect Covid-19 Test.**

**This product is for use under Emergency Use Authorization (EUA) only.**

**REF 21101  
MN-00005-US**

**Effective Date: October 2021  
Rev 5**

## Trademark, Patents and Copyright Statements

Detect™ is a trademark of Detect, Inc.

The purchase of this product conveys to the buyer the non-transferable right to use it in accordance with this user manual. No other rights are conveyed expressly, by implication or by estoppel. Furthermore, no rights for resale are conferred with the purchase of this product.

Copyright © Detect, Inc. 2021.  
All rights reserved.

This product has not been FDA cleared or approved, but has been authorized by FDA under an EUA.

This product has been authorized only for the detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens.

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.



Detect, Inc.  
530 Old Whitfield St.  
Guilford, CT 06437 USA  
+1 (855) 322-3692  
detect.com


# Contents

Using the Hub for the Detect Covid-19 Test	04
Product Information	06
Proprietary Name	06
Established Product Name	06
Intended Use	06
Hub Information	07
Hub Components	07
Device Diagram	07
Hub Specifications	08
Hub Installation	09
Hub Safety Precautions	10
Cleaning & Maintenance	11
Troubleshooting	12
Legend of Symbols	14
Repeat Use	15

# 1


## Using the Hub for the Detect Covid-19 Test

- 1.1 To start your test, you need a Detect Covid-19 Test (sold separately) and the Detect App (free to download). Follow the instructions in the app to prepare your sample using the Detect test. At the appropriate step, connect the Hub to a wall power outlet via the included USB-C cable and the USB power adapter. The Hub will beep when powered up, and it is ready to use when a solid green light is visible.

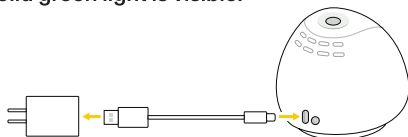

Caution: Do not insert the Test Tube into the Hub without powering up.

- 1.2 When ready, ensure that the Detect Cap on the Test Tube is securely tightened. Place the Test Tube in the Hub until you hear a beep. This indicates that the processing program has started.

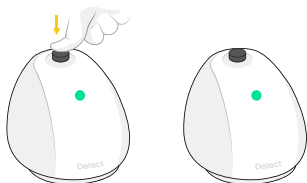

- 1.3 A blinking green LED light will indicate that the processing program is in progress.



- 1.4 A solid green LED light and two audible beeps will indicate completion of the process. This will take approximately 55 minutes.



- 1.5 Review the next steps in the Detect App before removing the Test Tube from the Hub.

## 2

## Product Information

- 2.1 Proprietary Name  
Detect Hub
- 2.2 Established Product Name  
Detect Hub
- 2.3 Intended Use  
This Hub is intended for use in combination with the Detect Covid-19 Test for the purpose of detecting SARS-CoV-2 nucleic acids in human anterior nares specimens.


## 3

## Hub Information

- 3.1 Hub Components  
Each Detect Hub comes with the following components:

<u>Component Name</u>	<u>Quantity</u>
Detect Hub .....	1
USB-C Cable .....	1
USB Power Adapter .....	1
Model S010WU0500200	
Input: 100-240V~50/60Hz 400mA	
Output: 5.0V=2000mA	
User Manual .....	1

- 3.2 Device Diagram

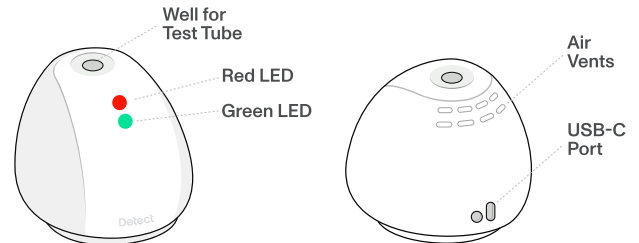


### 3.3 Hub Specifications

Model .....	21101
Input power .....	DC 5V==2A Max
Max. power .....	10 W
Operating environment .....	Indoor Use Only
Operating temperature range .....	59 °F - 86 °F (15 °C to 30 °C)
Storage temperature range .....	32 °F - 176 °F (0 °C to 80 °C)
Operating humidity range .....	25%–80%
Storage humidity range .....	0-80%
Max operating altitude .....	6500 ft (1981.2 m)
Dimensions (L×W×H) .....	4.3" x 3.4" x 2.6" (109 × 85 × 65 mm)
Net weight .....	3.3 oz (93 g)

### 3.4 Hub Installation

Only use the USB-C cable and USB power adapter included with this Detect Hub. Place the Hub on a level, horizontal surface in a dry place away from strong sunlight and with good ventilation around the device.

Attach the USB power adapter to the USB-C cable at the back of the Hub, and plug the other end of the power cord into an electrical outlet. The voltage requirement is 100-240V AC.

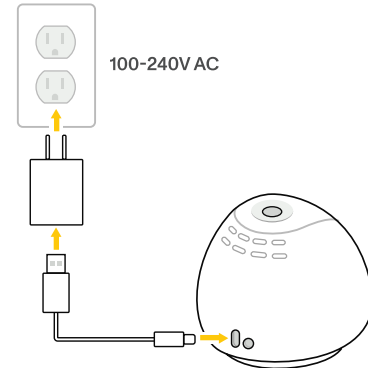



## 4

## Hub Safety Precautions



Please read this manual carefully before using the device and follow all recommended safety precautions to avoid personal injury or damage to the device:

- Only use with the included adapter.
- Do not place anything besides Detect Test Tubes into the Hub.
- Do not open the Test Tube during or after the processing step in order to prevent Hub contamination with biohazardous substance.
- Keep the Hub dry, dust-free and with good ventilation on all sides. Do not place near a heat source.
- Unplug the Detect Hub before cleaning.
- Do not attempt to open or repair the Hub. Contact Detect Customer Support for any questions.
- Power off the Hub when not in use by unplugging from the wall outlet.
- If it will be left unused for an extended period of time, keep the device covered with a cloth to protect it from dust.

## 5

## Cleaning & Maintenance



Power off the Hub before cleaning. Do not use corrosive cleaning liquids or fill the well with cleaning solutions.

Clean the plastic surface with cloth dampened with isopropyl alcohol. Do not clean the interior wall of the Hub well.

## 6 Troubleshooting

### PROBLEM

Hub does not power on  
(no LEDs are lit)



### SOLUTION

Check the following connections:

- Hub and USB-C cable
- USB-C cable and A/C adapter
- A/C adapter to outlet

*If problem persists, contact support.*

Red LED fades in/  
out continuously



Cause: Hub is not receiving correct power. *Find the original A/C adapter that came with the Hub package and use it.*

Both Red and Green LED fade in/  
out continuously



Cause: Hub detects that its internal and/or external operating temperature is outside the specified range.

*Relocate Hub to an environment where the temperature is 15 °C - 30 °C (59 °F - 86 °F). You may need to wait until Hub cools down to new temperature before it is operational again.*

If the problem persists, please contact support at:  
support@detect.com +1 (855) 322-3692

### PROBLEM

Red LED blinks continuously while Test Tube is inside Hub



### SOLUTION

Cause: Test Tube placed in the Hub prior to the Hub being plugged in. Take out the tube and re-insert. If power was interrupted during processing, discard the test as it is no longer valid—start a new test.

Red LED blinks 5 times after Test Tube removal

Cause: Test Tube removed before processing step is complete. The test is no longer valid—do not proceed.

Solid Red LED



Cause: Hub detects internal error. *Contact support.*

Red LED blinks very rapidly, accompanied by a continuous series of beeps

Cause: Hub detects a severe internal error. *Remove Hub from power immediately. Contact support.*

## Legend of Symbols



This symbol certifies that the electromagnetic interference from the device is under limits approved by the Federal Communications Commission.



This symbol indicates that the device has been tested to applicable standards by Underwriters Laboratory (UL).



This symbol indicates direct current (DC) voltage.



For In Vitro Diagnostics Use.



This symbol indicates that caution is necessary when operating the device in order to avoid undesirable consequences. See Section 4 of this manual to see the full list of safety precautions.



This symbol indicates the date on which the device was manufactured.



This symbol identifies the manufacturer of the device.



This symbol indicates that the instruction manual should be consulted.



This symbol indicates that the device should be kept in dry conditions.



This symbol indicates that the device is fragile and shall be handled with care.



This symbol indicates that the device must not be used if the package holding the device is damaged.



This symbol identifies the manufacturer's catalogue number.



This symbol identifies the manufacturer's serial number.



This symbol identifies that the electronics device meets the safety requirements specified for Class II equipment according to IEC 61140.

## Repeat Use

The Detect Hub is reusable and can be used with any additional Detect Tests.

However, if you need to dispose of it, please see the information below.

This electronics product contains substances which could damage the environment if not properly handled. As such, it is not possible to dispose of it in the regular trash.

Please check with your local government or electronics store to verify how to recycle your product. You can also check with your local electronics store to see if they will be able to recycle this product for you.

Depending on the solution chosen, fees may apply.

# Detect

Contact support at:  
[support@detect.com](mailto:support@detect.com)  
+1 (855) 322-3692