

SARS-CoV-2 Antigen (CoV2Ag)

For Use Under Emergency Use Authorization Only

For *in vitro* diagnostic use.

For prescription use only.

Assay for the Detection of SARS-CoV-2 Antigen

Current Revision and Date ^a	Rev. 05, 2022-05	
Product Name	ADVIA Centaur SARS-CoV-2 Antigen (CoV2Ag)	<div>REF 11207866</div> <div>(100 tests)</div> <div>REF 11207867</div> <div>(500 tests)</div>
Abbreviated Product Name	ADVIA Centaur CoV2Ag	
Test Name/ID	CoV2Ag	
Systems	ADVIA Centaur XP system ADVIA Centaur XPT system	
Materials Required but Not Provided	ADVIA Centaur CoV2Ag Calibrator (CoV2Ag CAL)	REF 11208047
	ADVIA Centaur CoV2Ag Quality Control (CoV2Ag QC)	REF 11207868
	Sample Inactivation Media	REF 11208651
	ADVIA Centaur Ancillary Probe Wash 1	REF 03395373
	ADVIA Centaur Probe Wash 3	REF 03333963
	ADVIA Centaur Wash 1 (2 x 1500 mL)	REF 01137199 (112351)
	ADVIA Centaur Wash 1 (2 x 2500 mL)	REF 03773025
Specimen Types	Anterior nasal swab	
Sample Volume	100 µL	
Measuring Interval	0.10–1000 Index	

^a A vertical bar in the page margin indicates technical content that differs from the previous version.

Intended Use

The ADVIA Centaur® SARS-CoV-2 Antigen (CoV2Ag) assay is a chemiluminescent immunoassay intended for the qualitative detection of SARS-CoV-2 nucleocapsid antigen in anterior nasal (AN) swab specimens collected in Siemens Healthineers Sample Inactivation Media, from individuals who are suspected of COVID-19 by their healthcare provider within the first 7 days of symptom onset, using the ADVIA Centaur® XP and ADVIA Centaur® XPT systems. Testing is limited to laboratories certified under the Clinical Laboratory Improvement Amendments of 1988 (CLIA), 42 U.S.C 263a, that meet the requirements to perform moderate or high complexity tests.

The ADVIA Centaur SARS-CoV-2 Antigen (CoV2Ag) assay does not differentiate between SARS-CoV and SARS-CoV-2.

Results are for the identification of SARS-CoV-2 nucleocapsid antigen. Antigen is generally detectable in anterior nasal swab specimens during the acute phase of infection. Positive results indicate the presence of viral antigens, but clinical correlation with patient history and other diagnostic information is necessary to determine infection status. Positive results do not rule out bacterial infection or co-infection with other viruses. The agent detected may not be the definite cause of disease. Laboratories within the United States and its territories are required to report all results to the appropriate public health authorities.

Negative results should be treated as presumptive and confirmation with a molecular assay, if necessary, for patient management may be performed. Negative results do not rule out SARS-CoV-2 infection and should not be used as the sole basis for treatment or patient management decisions, including infection control decisions. Negative results should be considered in the context of a patient's recent exposures, history, and the presence of clinical signs and symptoms consistent with COVID-19.

The ADVIA Centaur SARS-CoV-2 Antigen (CoV2Ag) assay is intended for use by trained clinical laboratory personnel specifically instructed and trained in *in vitro* diagnostic procedures and proper infection control procedures.

In the United States, the ADVIA Centaur SARS-CoV-2 Antigen (CoV2Ag) assay is only for use under the Food and Drug Administration's Emergency Use Authorization.

Summary and Explanation

COVID-19 (coronavirus disease 2019) is the illness resulting from infection with SARS-CoV-2 (severe acute respiratory syndrome coronavirus 2) virus.¹⁻⁵ The virus spreads readily from person to person primarily through infected secretions, such as saliva and respiratory droplets or aerosols.⁶ Evidence supports spread by both symptomatic and asymptomatic individuals.⁷ The virus incubation period ranges from 2–14 days following exposure, with most cases showing symptoms within approximately 5 days after exposure.⁸

SARS-CoV-2 nucleic acid amplification testing, such as reverse transcription polymerase chain reaction (RT-PCR), is considered the gold standard for diagnostic testing for current infection, typically using an upper respiratory tract specimen. RT-PCR detects the genetic material of the virus, while antigen tests detect a viral protein (e.g. nucleocapsid). Immunoassays that detect the SARS-CoV-2 nucleocapsid antigen are also used for the diagnosis of current infection.⁹ Due to the highly contagious nature of this virus and the mode of transmission, the need for additional, less complex testing solutions has been recognized. Antigen testing, as part of a SARS-CoV-2 testing strategy, can be used to identify symptomatic and asymptomatic individuals currently infected with SARS-CoV-2.⁹

Principles of the Procedure

The ADVIA Centaur CoV2Ag assay is a fully automated 1-step sandwich immunoassay using acridinium ester chemiluminescent technology. The Solid Phase reagent containing streptavidin magnetic latex particles coated with biotinylated anti-SARS-CoV-2 mouse monoclonal capture antibodies and the Lite Reagent containing acridinium-ester-labeled anti-SARS-CoV-2 mouse monoclonal antibodies are added to the sample. The SARS-CoV-2 antigen in the sample complexes with the antibodies. Antibody-antigen complexes will form if SARS-CoV-2 antigen is present in the sample.

A direct relationship exists between the amount of SARS-CoV-2 nucleocapsid antigen present in the patient sample and the amount of relative light units (RLUs) detected by the system.

A result of reactive or nonreactive is determined according to the Index Value established with the calibrators. Refer to *Interpretation of Results*.

Reagents

Material Description	Storage	Stability
ADVIA Centaur CoV2Ag ReadyPack® primary reagent pack^{a, b} Lite Reagent 4.0 mL/reagent pack Mouse monoclonal anti-SARS-CoV-2 antibodies labeled with acridinium ester (~0.60 µg/mL); buffer; surfactant; bovine serum albumin (BSA); goat serum; mouse IgG; bovine gamma globulin; goat gamma globulin; sodium azide (< 0.1%) Solid Phase 10.5 mL/reagent pack Streptavidin-coated paramagnetic microparticles preformed with biotinylated SARS-CoV-2 monoclonal antibodies (~0.016 mg/mL); buffer; BSA; goat serum; bovine gamma globulin; surfactant; sodium azide (< 0.1%)	Unopened at 2–8°C Onboard	Until expiration date on product 21 days
ADVIA Centaur CoV2Ag ReadyPack ancillary reagent pack^{a, b} Ancillary Reagent 10.0 mL/reagent pack Buffer; BSA; preservatives	Unopened at 2–8°C Onboard	Until expiration date on product 21 days
ADVIA Centaur Ancillary Probe Wash 1 ReadyPack ancillary reagent pack^{a, c} 25.0 mL/pack 0.4 N sodium hydroxide	Unopened at 2–8°C Onboard	Until expiration date on product 14 days
ADVIA Centaur Probe Wash 3^{a, c} 50.0 mL/pack Sodium hypochlorite (0.5%); sodium hydroxide (< 0.5%); pH 11.0	Unopened at 2–8°C Onboard	Until expiration date on product 100 days
ADVIA Centaur Wash 1^{a, c} 1500 mL/pack Phosphate-buffered saline; sodium azide (< 0.1%); surfactant	Unopened at 2–25°C Onboard	Until expiration date on product 1 month

Material Description	Storage	Stability
ADVIA Centaur Wash 1 ^{a, c} 2500 mL/pack Phosphate-buffered saline; sodium azide (< 0.1%); surfactant	Unopened at 2–25°C Onboard	Until expiration date on product 1 month

^a Store in an upright position.

^b Prevent exposure to sunlight and heat.

^c Reagent is required, but not provided in the assay kit. Refer to *Materials Required but Not Provided*.

Warnings and Precautions

For Use Under Emergency Use Authorization Only

For *in vitro* diagnostic use only.

For prescription use only.

In the U.S., this product has not been FDA cleared or approved, but has been authorized by FDA under an Emergency Use Authorization (EUA) for use by authorized laboratories; use by laboratories certified under the Clinical Laboratory Improvement Amendments of 1988 (CLIA), 42 U.S.C 263a, that meet the requirements to perform moderate or high complexity tests.

This product has been authorized only for the detection of proteins from SARS-CoV-2, not for any other viruses or pathogens.

In the U.S., the emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of *in vitro* diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug, and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or the authorization is revoked sooner.

Laboratories within the U.S. and its territories are required to report all results to the appropriate public health laboratories.

Do not use this kit beyond the use by date printed on the product label.

Do not pool the contents of different vials of the same reagent (even if the reagents are from the same lot). Calibrators are also lot-specific and cannot be shared between kits.

Treat all specimens as potentially infectious. Follow universal precautions when handling samples, this kit, and its contents.

Proper laboratory safety techniques should be followed at all times when working with SARS-CoV-2 patient samples.

All results must be interpreted together with other clinical information available.

In order to obtain accurate results, the test must follow this package insert.

Inadequate or inappropriate sample collection, storage, and transport can result in incorrect results.

CAUTION

Federal (USA) law restricts this device to sale by or on the order of a licensed healthcare professional.

Safety data sheets (SDS) available on [siemens-healthineers.com](https://www.siemens-healthineers.com).



H290, H319, H315
P264, P280,
P305+P351+P338

Warning!

May be corrosive to metals. Causes serious eye irritation. Causes skin irritation.

Wash hands thoroughly after handling. Wear protective gloves/protective clothing/eye protection/face protection. IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.

Contains: sodium hydroxide (ADVIA Centaur Ancillary Probe Wash 1)

CAUTION

This device contains material of animal origin and should be handled as a potential carrier and transmitter of disease.

Contains sodium azide as a preservative. Sodium azide can react with copper or lead plumbing to form explosive metal azides. On disposal, flush reagents with a large volume of water to prevent buildup of azides. Disposal into drain systems must be in compliance with prevailing regulatory requirements.

Do not eat, drink, or smoke in the area where the specimens and kit contents are handled. Use appropriate precautions in the collection, handling, storage, and disposal of patient samples and used kit contents.

Dispose of hazardous or biologically contaminated materials in accordance with federal, state, and local requirements. Discard all materials in a safe and acceptable manner and in compliance with prevailing regulatory requirements.

Storage and Stability

Store all reagents at 2–8°C in an upright position, away from light and heat. Do not use products beyond the expiration date printed on the product labeling. Do not use if the test device package is damaged.

For information about product storage and stability, refer to *Reagents*.

Onboard Stability

Discard products at the end of the onboard stability interval. Do not use products beyond the expiration date printed on the product labeling.

For information about product onboard stability, refer to *Reagents*.

Specimen Collection and Handling

Anterior nasal (AN) swab specimens stored in Sample Inactivation Media ([REF](#) 11208651) is the recommended sample type for this assay. This assay has not been authorized for use with transport media other than Sample Inactivation Media.

Inadequate or inappropriate sample collection, storage, and transport can result in incorrect results.

Do not use heat-inactivated specimens.

Use appropriate precautions when collecting, handling, storing, and transporting swab specimens. Refer to the CDC's "Interim Guidelines for Collection, Handling and Transportation of clinical specimens from persons with Coronavirus Disease 2019 (COVID-19)"¹⁰ and to the WHO's interim guidance document, "Laboratory testing for coronavirus disease (COVID-19) in suspected human cases"¹¹, for proper procedures.

The handling and storage information provided here is based on data or references maintained by the manufacturer. It is the responsibility of the individual laboratory to use all available references and/or its own studies when establishing alternate stability criteria to meet specific needs.

Collecting the Specimen

- Handle all specimens as if they are capable of transmitting disease.¹²
- When collecting any type of sample, follow universal collection precautions and guidelines according to your organization. Users should be trained in appropriate sample collection and handling procedures.
- For the collection of AN swab specimens, follow the Centers for Disease Control and Prevention (CDC) Swab Collection Guidelines¹⁰ and the swab manufacturers' recommendations.
- Keep tubes capped at all times.

Storing the Specimen

- Inactivated samples are stable for up to 4 hours onboard the system.
- Inactivated samples can be stored at 2–8°C for up to 3 days.
- Freeze inactivated samples at $\leq -20^{\circ}\text{C}$ for longer storage. Do not store in a frost-free freezer. Avoid more than 2 freeze-thaw cycles. Thoroughly mix thawed inactivated samples and centrifuge them before using.

Transporting the Specimen

Package and label specimens for shipment in compliance with applicable federal and international regulations covering the transport of clinical specimens and etiological agents.

Ship specimens frozen if specimen storage exceeds sample stability specified in *Storing the Specimen*.

Preparing the Samples

This assay requires 100 μL of sample for a single determination. This volume does not include the unusable volume in the sample container or the additional volume required when performing duplicates or other tests on the same sample. For a complete list of appropriate sample containers and information about determining the minimum required volume, refer to the system online help.

Do not use samples with apparent contamination.

Sample Preparation of Swab Specimens Prior to Testing

Perform the following steps to prepare specimens received in Sample Inactivation Media (REF 11208651) prior to loading the sample on the system.

WARNING

Use appropriate precautions for handling potentially infectious specimens when performing the procedures described below. Refer to *Specimen Collection and Handling*.

-
1. Allow the sample to thaw, if necessary.
 2. Carefully remove the cap from the sample tube. Roll the swab at least 5 times against the side of the tube and remove the swab. Dispose of the used swab into a biohazardous waste collection container.

- The inactivated sample should be tested as soon as possible. If the sample cannot be tested immediately, store the sample according to the storage guidance in *Storing the Specimen*.

Note Transferring the inactivated sample received in Sample Inactivation Media (REF 11208651) to another cup is not required, as it can be placed directly on the ADVIA Centaur XP and ADVIA Centaur XPT systems.

- Before placing samples on the system, ensure that samples are free of bubbles, foam, or particulate matter. Remove particulates by centrifugation, if necessary.
- Place the sample cup containing the inactivated sample on the system for testing.

Procedure

Materials Provided

The following materials are provided:

REF	Contents	Number of Tests
11207866	1 ReadyPack primary reagent pack containing ADVIA Centaur CoV2Ag Lite Reagent and Solid Phase 1 ReadyPack ancillary reagent pack containing ADVIA Centaur CoV2Ag Ancillary Reagent (ANC) ADVIA Centaur CoV2Ag master curve card	100
11207867	5 ReadyPack primary reagent packs containing ADVIA Centaur CoV2Ag Lite Reagent and Solid Phase 5 ReadyPack ancillary reagent packs containing ADVIA Centaur CoV2Ag Ancillary Reagent (ANC) ADVIA Centaur CoV2Ag master curve card	500

Materials Required but Not Provided

The following materials are required to perform these assays, but are not provided:

REF	Description	
	ADVIA Centaur XP system ^a ADVIA Centaur XPT system ^a	
11208047	ADVIA Centaur CoV2Ag CAL (calibrator)	2 x 1.0 mL low calibrator (CAL L) 2 x 1.0 mL high calibrator (CAL H) calibrator assigned value sheet
11207868	ADVIA Centaur CoV2Ag QC (quality control material)	2 x 1.6 mL negative quality control, level 1 (CONTROL - 1) 2 x 1.6 mL positive quality control, level 2 (CONTROL + 2) quality control assigned value sheet
11208651	Sample Inactivation Media (inactivation reagent)	50 x 1.0 mL
03395373	ADVIA Centaur Ancillary Probe Wash 1 (probe wash)	2 ReadyPack ancillary reagent packs containing 25.0 mL/ pack (APW 1)
03333963	ADVIA Centaur Probe Wash 3 (probe wash)	50.0 mL (PW 3)

REF	Description	
01137199 (112351)	ADVIA Centaur Wash 1 (wash)	2 x 1500 mL/pack WASH 1
03773025	ADVIA Centaur Wash 1 (wash)	2 x 2500 mL/pack WASH 1

^a Additional system fluids are required to operate the system: ADVIA Centaur Acid Reagent, ADVIA Centaur Base Reagent, and ADVIA Centaur Cleaning Solution.

Assay Procedure

The system automatically performs the following steps:

1. Dispenses 100 µL of sample into a cuvette.
2. Dispenses 60 µL of Ancillary Reagent, then incubates for 6 minutes at 37°C.
3. Dispenses 105 µL of Solid Phase and 40 µL of Lite Reagent, then incubates for 18 minutes at 37°C.
4. Separates the Solid Phase from the mixture and aspirates the unbound reagent.
5. Performs a wash sequence using ADVIA Centaur Wash 1.
6. Dispenses 300 µL each of ADVIA Centaur Acid Reagent and ADVIA Centaur Base Reagent to initiate the chemiluminescent reaction.
7. Reports results.

Preparing the Reagents

All reagents are liquid and ready to use. Before loading the packs onto the system, reagents require mixing. For information about mixing the reagents, refer to the system online help.

Note The Ancillary Reagent provided in this kit is matched to the Solid Phase and Lite Reagent. Do not mix Ancillary Reagent lots with different lots of Solid Phase and Lite Reagent.

Preparing the System

Ensure that sufficient materials are loaded on the system. Refer to *Materials Provided* and *Materials Required but Not Provided* for guidance about required reagents.

For information about loading products, refer to the system online help.

Master Curve Definition

Before initiating calibration on each new lot of reagent, enter the assay master curve values by scanning the master curve card. For information about defining the master curve, refer to the system online help.

Performing Calibration

For calibration of the ADVIA Centaur CoV2Ag assay, use the ADVIA Centaur CoV2Ag CAL. Use the calibrators in accordance with the calibrator instructions for use.

Calibration Frequency

Perform a calibration if one or more of the following conditions exist:

- At the end of the 21-day calibration interval.
- When changing lot numbers of primary reagent packs.
- When indicated by quality control results.
- After major maintenance or service, if indicated by quality control results.

Follow government regulations or accreditation requirements for calibration frequency. Individual laboratory quality control programs and procedures may require more frequent calibration.

Performing Quality Control

For quality control of the ADVIA Centaur CoV2Ag assay, use the ADVIA Centaur CoV2Ag QC or an equivalent product at least once during each day that samples are analyzed. Use the quality control material in accordance with the quality control instructions for use. For the assigned values, refer to the quality control assigned value sheet provided.

In addition, perform quality control:

- Following a valid calibration
- With use of a new lot of reagent
- When troubleshooting test results that do not match clinical conditions or symptoms

Follow government regulations or accreditation requirements for quality control frequency. Individual laboratory quality control programs and procedures may require more frequent quality control testing.

Acceptable performance is achieved when the analyte values obtained are within the expected control interval for the system, as indicated by the manufacturer of the control material or within the interval determined by an internal laboratory quality control procedure.

Follow your laboratory's quality control procedures if the results obtained do not fall within the acceptable limits. For information about entering quality control definitions, refer to the system online help.

Taking Corrective Action

If the quality control results do not fall within the expected control interval, do not report results. Perform corrective actions in accordance with established laboratory protocol. For suggested protocol, refer to the system online help.

Results

Calculation of Results

The system determines the result using the calculation procedure described in the system online help. The system reports results in Index Values. Refer to *Interpretation of Results*.

Interpretation of Results

The system reports ADVIA Centaur CoV2Ag assay results in Index Values and as Nonreactive or Reactive:

- **Nonreactive:** < 1.00 Index. These samples are considered presumptive negative for SARS-CoV-2 nucleocapsid antigen.
- **Reactive:** ≥ 1.00 Index. These samples are considered positive for SARS-CoV-2 nucleocapsid antigen.

The cut-off value for the ADVIA Centaur CoV2Ag assay was verified based on clinical agreement of results.

Results of this assay should always be interpreted in conjunction with the patient's medical history, clinical presentation, and other findings.

Limitations

The following information pertains to limitations of the assay:

- This device has been evaluated for use with human specimen material only.
- The performance of the assay has not been established with specimen types other than those defined in the *Intended Use*.
- The performance of this test was evaluated using Sample Inactivation Media only. The use of other viral transport media should not be used in conjunction with this test.
- The performance of this test has not yet been clinically validated for use in patients without signs and symptoms of respiratory infection.
- Failure to follow the instructions for use may adversely affect test performance and/or invalidate the test result.
- Results should be considered in the context of a patient's recent exposures, history, and the presence of clinical signs and symptoms consistent with COVID-19, and confirmed with a molecular assay, if necessary, for patient management.
- A false negative result may occur due to inadequate sample collection, storage, and/or handling.
- A false negative result may occur if the concentration of the target antigen in the clinical specimen is below the detection limits of the device.
- Positive test results do not rule out co-infections with other pathogens.
- Negative results should be treated as presumptive and confirmation with a molecular assay, if necessary for patient management, may be performed.
- Percent interference > 15% was observed in contrived samples of 100% pooled human nasal wash, which may result in a false negative.
- Percent interference > 15% was observed in contrived samples containing 15% Afrin, which may result in a false negative.
- The performance of this test is dependent upon proper specimen collection, handling, storage, and preparation. Failure to observe proper procedures in any one of these steps can lead to incorrect results.
- Results from antigen testing should not be used as the sole basis to diagnose or exclude SARS-CoV-2 infection or to determine infection status.
- This test does not differentiate between SARS-CoV and SARS-CoV-2.
- If the differentiation of specific SARS viruses and strains is needed, additional testing, in consultation with state or local public health departments, is required.
- This test cannot rule out diseases caused by other bacterial or viral pathogens.
- This test will indicate the presence of SARS-CoV-2 nucleocapsid protein antigen in the specimen from both viable and non-viable SARS-CoV-2 virus. Test performance depends on the amount of virus (antigen) in the sample and may or may not correlate with viral culture results performed on the same sample.
- This device is a qualitative test and does not provide information on the viral concentration present in the specimen.

- The performance of this test was established based on the evaluation of a limited number of clinical specimens collected between January 25, 2022 and February 10, 2022. The clinical performance has not been established in all circulating variants but is anticipated to be reflective of the prevalent variants in circulation at the time and location of the clinical evaluation. Performance at the time of testing may vary depending on the variants circulating, including newly emerging strains of SARS-CoV-2 and their prevalence, which change over time.
- Positive and negative predictive values are highly dependent on prevalence. False-negative test results are more likely during peak activity when the prevalence of the disease is high. False-positive test results are more likely during the periods of low SARS-CoV-2 activity when prevalence is moderate to low.

Conditions of Authorization for the Laboratory

The ADVIA Centaur SARS-CoV-2 Antigen (CoV2Ag) assay Letter of Authorization, along with the authorized Fact Sheet for Healthcare Providers, the authorized Fact Sheet for Patients, and authorized labeling are available on the FDA website:

<https://www.fda.gov/medical-devices/coronavirus-disease-2019-covid-19-emergency-use-authorizations-medical-devices/in-vitro-diagnostics-euas>

Authorized laboratories using the ADVIA Centaur CoV2Ag assay must adhere to the Conditions of Authorization indicated in the Letter of Authorization as listed below:

- Authorized laboratories^a using this product must include with test result reports, all authorized Fact Sheets. Under exigent circumstances, other appropriate methods for disseminating these Fact Sheets may be used, which may include mass media.
- Authorized laboratories using this product must use the product as outlined in the authorized labeling. Deviations from the authorized procedures, including the authorized instruments, authorized clinical specimen types, authorized control materials, authorized other ancillary reagents and authorized materials required to use the ADVIA Centaur CoV2Ag assay are not permitted.
- Authorized laboratories that receive this product must notify the relevant public health authorities of their intent to run the assay prior to initiating testing.
- Authorized laboratories using this product must have a process in place for reporting test results to healthcare providers and relevant public health authorities, as appropriate.
- Authorized laboratories must collect information on the performance of this product and report to DMD/OHT7-OIR/OPEQ/ CDRH (via email: CDRH_EUA_Reporting@fda.hhs.gov) and Siemens Healthineers Technical Support (<https://www.siemens-healthineers.com/en-us/>; tel: 1-877-229-3711) any suspected occurrence of false reactive (positive) or false nonreactive (negative) results and significant deviations from the established performance characteristics of the assay of which they become aware.

- All laboratory personnel using this product must be appropriately trained in automated immunoassay techniques and use appropriate laboratory and personal protective equipment when handling this kit, and use this product in accordance with the authorized labeling. All laboratory personnel using the assay must also be trained in and be familiar with the interpretation of results of this product.
- Siemens Healthineers, authorized distributors, and authorized laboratories using this product must ensure that any records associated with this EUA are maintained until otherwise notified by FDA. Such records will be made available to FDA for inspection upon request.

^a The letter of authorization refers to, "Laboratories certified under the Clinical Laboratory Improvement Amendments of 1988 (CLIA), 42 U.S.C. §263a, that meet the requirements to perform moderate or high complexity tests" as "authorized laboratories".

Performance Characteristics

The reagent formulations used on the ADVIA Centaur XP and ADVIA Centaur XPT systems are the same as those used on the Atellica® IM Analyzer. Some performance characteristics for the ADVIA Centaur CoV2Ag assay were established with the Atellica IM CoV2Ag assay using the Atellica IM Analyzer.

Assay performance characteristics are representative data. Results obtained at individual laboratories may vary from the data presented.

Measuring Interval

0.10–1000 Index is reported as Nonreactive (< 1.00 Index) or Reactive (≥ 1.00 Index).

Detection Capability

Limit of Detection (LoD) = 31.2 TCID₅₀/mL

LOD was determined by evaluating different dilutions of gamma irradiated viral lysate SARS-CoV-2 (USA WA1/2020) in Sample Inactivation Media ([REF](#) 11208651). The study was performed by adding 50µL of viral solution to dry swabs and eluting the swabs in Sample Inactivation Media ([REF](#) 11208651). LoD is defined as the lowest virus concentration at which a minimum of 19 out of 20 generate a Reactive result. Testing was performed with the ADVIA Centaur CoV2Ag assay using the ADVIA Centaur XP system.

The LoD corresponds to the lowest concentration of SARS-CoV-2 antigen that can be detected with a probability of 95%.

Clinical Agreement

Positive percent agreement and negative percent agreement were determined by testing a total of 165 anterior nasal (AN) swab samples in accordance with CLSI Document EP12-A2.¹³ Samples were collected from 68 male and 97 female symptomatic individuals suspected of COVID-19 and enrolled sequentially (all-comers), with one or more confirmed symptoms based on CDC guidelines and within 7 days of symptom onset. All samples were collected between January and February 2022 from 2 vendors at 5 sites in the U.S. Of the 47 RT-PCR positive samples, 44 were confirmed as the Omicron variant using next-generation sequencing and 3 samples failed sequencing due to an insufficient quantity of viral material.

The samples were collected in Sample Inactivation Media ([REF](#) 11208651), refrigerated from the time of collection, and tested with the ADVIA Centaur CoV2Ag assay using the ADVIA Centaur XP system. The results were compared to FDA EUA-authorized RT-PCR comparator results.

Subject Demographics by Age

Age (Years)	Anterior Nasal Swab Samples (N=165)		
	Number Tested	ADVIA Centaur CoV2Ag Assay - Reactive	Prevalence
< 22	29	9	31%
22–59	117	29	25%
≥ 60	19	2	11%

Performance in Anterior Nasal Swab Samples Collected from RT-PCR Positive Symptomatic and RT-PCR Negative Individuals

		Comparative PCR Method		
		Positive	Negative	Total
ADVIA Centaur CoV2Ag Assay	Reactive	40	0	40
	Nonreactive	7	118	125
	Total	47	118	165
Positive Percent Agreement: 85.11% (95% Confidence Interval: 72.31%–95.59%)				
Negative Percent Agreement: 100% (95% Confidence Interval: 96.85%–100%)				

Performance in RT-PCR Positive Anterior Nasal Swab Samples Collected from Symptomatic Individuals from Date of Symptom Onset

Days Since Symptom Onset	Comparative PCR Method - Positive (Cumulative)	ADVIA Centaur CoV2Ag Assay - Reactive (Cumulative)	Positive Percent Agreement (%)	95% Confidence Interval
1	11	8	72.73%	43.44%–90.25%
2	22	18	81.82%	61.48%–92.69%
3	30	24	80.00%	62.69%–90.49%
4	35	28	80.00%	64.11%–89.96%
5	40	33	82.50%	68.05%–91.25%
6	44	37	84.09%	70.63%–92.07%
7	47	40	85.11%	72.31%–92.59%

Interferences

Interference testing was performed using the Atellica IM CoV2Ag assay with the Atellica IM Analyzer. Percent interference > 15% was observed for Afrin (oxymetazoline) and pooled human nasal wash; this did not result in a false negative result. None of the other substances tested were found to interfere with the clinical interpretation of the test in nonreactive and weakly reactive samples at the substance test concentrations indicated. The following results are representative of the performance of the assay:

Substance	Substance Test Concentration	Nonreactive Sample	Spiked Reactive Sample
Afrin (oxymetazoline)	15% v/v	Negative	Positive
Alkalol homeopathic nasal wash	1:10 dilution	Negative	Positive
Biotin	3500 ng/mL	Negative	Positive
Chloraseptic (benzocaine/menthol)	1.5 mg/mL	Negative	Positive
CVS Nasal Drops (phenylephrine)	15% v/v	Negative	Positive
CVS Nasal Spray (cromolyn)	15% v/v	Negative	Positive
Fluticasone propionate	5% v/v	Negative	Positive
Hemoglobin	75 mg/dL	Negative	Positive
Mucin	0.5%	Negative	Positive
Mupirocin	10 mg/mL	Negative	Positive
Naso GEL (NeilMed)	5% v/v	Negative	Positive
Pooled human nasal wash	100% v/v	Negative	Positive
Sore throat phenol spray	15% v/v	Negative	Positive
Tamiflu (oseltamivir phosphate)	5 mg/mL	Negative	Positive
Tobramycin	4 µg/mL	Negative	Positive
Whole blood	4% v/v	Negative	Positive
Zicam Nasal Spray	5% v/v	Negative	Positive

Cross-Reactivity

The assay was evaluated for potential cross-reactivity in specimens with various microorganisms and viruses using the Atellica IM CoV2Ag assay with the Atellica IM Analyzer. The following results are representative of the performance of the assay:

Microorganism/Virus	Concentration Tested	Nonreactive Sample	Spiked Reactive Sample
Adenovirus 7a	10 ⁵ PFU/mL ^a	Negative	Positive
<i>Bordetella pertussis</i>	10 ⁶ CFU/mL ^b	Negative	Positive
<i>Candida albicans</i>	10 ⁶ CFU/mL	Negative	Positive
<i>Chlamydia pneumoniae</i>	10 ⁶ IFU/mL ^c	Negative	Positive
Enterovirus 68	10 ⁵ PFU/mL	Negative	Positive

Microorganism/Virus	Concentration Tested	Nonreactive Sample	Spiked Reactive Sample
<i>Haemophilus influenzae</i>	10 ⁶ CFU/mL	Negative	Positive
Human coronavirus 229E	10 ⁵ PFU/mL	Negative	Positive
Human coronavirus NL63	10 ⁵ PFU/mL	Negative	Positive
Human coronavirus OC43	10 ⁵ PFU/mL	Negative	Positive
Human Metapneumovirus (hMPV)	10 ⁵ PFU/mL	Negative	Positive
Human parainfluenza virus 1 (HPIV-1)	10 ⁵ PFU/mL	Negative	Positive
Human parainfluenza virus 2 (HPIV-2)	10 ⁵ PFU/mL	Negative	Positive
Human parainfluenza virus 3 (HPIV-3)	10 ⁵ PFU/mL	Negative	Positive
Human parainfluenza virus 4 (HPIV-4)	10 ⁵ PFU/mL	Negative	Positive
Influenza A (H1N1)	10 ⁵ PFU/mL	Negative	Positive
Influenza B	10 ⁵ PFU/mL	Negative	Positive
<i>Legionella pneumophila</i>	10 ⁶ CFU/mL	Negative	Positive
Middle East respiratory syndrome coronavirus (MERS-CoV)	10 ⁵ PFU/mL	Negative	Positive
<i>Mycobacterium tuberculosis</i>	10 ⁶ CFU/mL	Negative	Positive
<i>Mycoplasma pneumoniae</i>	10 ⁶ CCU/mL ^d	Negative	Positive
<i>Pneumocystis jirovecii</i> (PJP)	10 ⁶ CFU/mL	Negative	Positive
Respiratory syncytial virus (RSV)	10 ⁵ PFU/mL	Negative	Positive
Rhinovirus	10 ⁵ PFU/mL	Negative	Positive
<i>Staphylococcus aureus</i>	10 ⁶ CFU/mL	Negative	Positive
<i>Staphylococcus epidermidis</i>	10 ⁶ CFU/mL	Negative	Positive
<i>Streptococcus pneumoniae</i>	10 ⁶ CFU/mL	Negative	Positive
<i>Streptococcus pyogenes</i>	10 ⁶ CFU/mL	Negative	Positive
<i>Streptococcus salivarius</i>	10 ⁶ CFU/mL	Negative	Positive

^a PFU = plaque-forming unit

^b CFU = colony-forming unit

^c IFU = inclusion forming unit

^d CCU = color changing unit

In Silico analysis using the Basic Local Alignment Search Tool (BLAST) managed by the National Center for Biotechnology Information (NCBI) was used to assess the degree of protein sequence homology to estimate the cross-reactivity with SARS-CoV-2 virus in the presence of organisms that were not available for wet testing. There was 36.74% homology between SARS-CoV-2 nucleocapsid and human coronavirus HKU1 nucleocapsid, suggesting a low likelihood of cross-reactivity; however, cross-reactivity cannot be ruled out. The comparison between SARS-CoV-2 nucleocapsid protein and SARS-CoV-1 shows homology of 90.52% and suggests that there may be significant cross-reactivity in this assay.

High-Dose Hook Effect

No false negative results due to a high-dose hook effect were found with the ADVIA Centaur CoV2Ag assay up to 2.6×10^5 TCID₅₀/mL.

Standardization

The assay standardization for the ADVIA Centaur CoV2Ag assay is based on agreement with known positive and negative SARS-CoV-2 samples. Assigned values for calibrators are traceable to this standardization.

Technical Assistance

For customer support, contact your local technical support provider or distributor.

siemens-healthineers.com

References






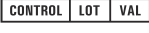
1. Lai CC, Shih TP, Ko WC, Tang HJ, Hsueh PR. Severe acute respiratory syndrome coronavirus 2 (SARS-CoV-2) and coronavirus disease-2019 (COVID-19): The epidemic and the challenges. *Int J Antimicrob Agents*. 2020;55(3):105924.
2. Guo YR, Cao QD, Hong ZS, Tan YY, Chen SD, Jin HJ, *et al*. The origin, transmission and clinical therapies on coronavirus disease 2019 (COVID-19) outbreak - an update on the status. *Mil Med Res*. 2020;7(1):11.
3. Li H, Liu SM, Yu XH, Tang SL, Tang CK. Coronavirus disease 2019 (COVID-19): current status and future perspectives. *Int J Antimicrob Agents*. 2020;29:105951.
4. Wu Z and McGoogan JM. Characteristics of and Important Lessons From the Coronavirus Disease 2019 (COVID-19) Outbreak in China: Summary of a Report of 72314 Cases From the Chinese Center for Disease Control and Prevention. *JAMA* 2020;323(13):2648.
5. Pan A, Liu L, Wang C, Guo H, Hao X, Wang Q, *et al*. Association of Public Health Interventions With the Epidemiology of the COVID-19 Outbreak in Wuhan, China. *JAMA*. 2020:6130.
6. World Health Organization. Transmission of SARS-CoV-2: implications for infection prevention precautions. Scientific Brief. WHO Reference Number: WHO/2019-nCoV/Sci_Brief/Transmission_modes/2020.3. Published July 9, 2020. Accessed January 13, 2021.
7. Lai CC, Liu YH, Wang CY, Wang YH, Hsueh SC, Yen MY, *et al*. Asymptomatic carrier state, acute respiratory disease, and pneumonia due to severe acute respiratory syndrome coronavirus 2 (SARS-CoV-2): Facts and myths. *J Microbiol Immunol Infect*. 2020;pii: S1684-1182;(20):30040-2.
8. Lauer S, Grantz K, Bi Q, *et al*. The Incubation Period of Coronavirus Disease 2019 (COVID-19) From Publicly Reported Confirmed Cases: Estimation and Application. *Annals of Internal Medicine*. 2020 May 5;172(9):577-582. doi:10.7326/M20-0504.
9. Centers for Disease Control and Prevention. Interim Guidance for Antigen Testing for SARS-CoV-2. <https://www.cdc.gov/coronavirus/2019-ncov/lab/resources/antigen-tests-guidelines.html>. Accessed Jan. 13, 2021.
10. Centers for Disease Control and Prevention. Interim Guidelines for Collecting, Handling, and Testing Clinical Specimens for COVID-19. <https://www.cdc.gov/coronavirus/2019-ncov/lab/guidelines-clinical-specimens.html>. Accessed Jan. 13, 2021.
11. World Health Organization. Laboratory testing for coronavirus disease (COVID-19) in suspected human cases. <https://www.who.int/publications/i/item/laboratory-testing-for-2019-novel-coronavirus-in-suspected-human-cases-20200117>. Published March 19, 2020. Accessed January 13, 2021.

12. Clinical and Laboratory Standards Institute. *Protection of Laboratory Workers From Occupationally Acquired Infections; Approved Guideline—Fourth Edition*. Wayne, PA: Clinical and Laboratory Standards Institute; 2014. CLSI Document M29-A4.
13. Clinical and Laboratory Standards Institute. *User Protocol for Evaluation of Qualitative Test Performance; Approved Guideline—Second Edition*. Wayne, PA: Clinical and Laboratory Standards Institute; 2008. CLSI Document EP12-A2.

Definition of Symbols

The following symbols may appear on the product labeling:

Symbol	Symbol Title	Symbol	Symbol Title
	Manufacturer		Authorized representative in the European Community
	Use-by date		Batch code
	Catalog number		Contains sufficient for <n> tests
	Consult Instructions for Use		Version of Instructions for Use
	Internet URL address to access the electronic instructions for use		Revision
	In vitro diagnostic medical device		Unique Device Identifier
	Prescription device (US only)		CE Marking
	CE Marking with Notified Body		Keep away from sunlight
	Temperature limit		Lower limit of temperature
	Upper limit of temperature		Do not freeze
	Do not re-use		This way up
	Recycle		Caution
	Biological risks		Document face up ^a
	Common Units		International System of Units
	Date format (year-month-day)		Date format (year-month)
	Target		Interval

Symbol	Symbol Title	Symbol	Symbol Title
	Handheld barcode scanner		Variable hexadecimal number that ensures the Master Curve and Calibrator definition values entered are valid.
	Lot details		Master Curve definition
	Calibrator lot value		Quality control lot value

^a Indicates Assay-eNote

Legal Information

Atellica, ReadyPack, and ADVIA Centaur are trademarks of Siemens Healthineers.

© 2022 Siemens Healthineers. All rights reserved.



Siemens Healthcare Diagnostics Inc.
511 Benedict Avenue
Tarrytown, NY 10591 USA

Siemens Healthineers Headquarters

Siemens Healthcare GmbH
Henkestraße 127
91052 Erlangen
Germany
Phone: +49 9131 84-0
siemens-healthineers.com

SARS-CoV-2 Antigen Calibrator (CoV2Ag CAL)

Current Revision and Date^a	Rev. 03, 2022-05
Product Name	ADVIA Centaur SARS-CoV-2 Antigen Calibrator (CoV2Ag CAL)
Abbreviated Product Name	ADVIA Centaur CoV2Ag CAL
	2 x 1.0 mL low calibrator CAL L REF 11208047 2 x 1.0 mL high calibrator CAL H Calibrator assigned value sheets and barcode labels CAL LOT VAL
Systems	ADVIA Centaur systems

^a A vertical bar in the page margin indicates technical content that differs from the previous version.

FOR USA:

For Use Under Emergency Use Authorization Only

For *in vitro* diagnostic use.

For Prescription Use.

For Professional Use.

Intended Use

The ADVIA Centaur® SARS-CoV-2 Antigen Calibrator (CoV2Ag CAL) is for *in vitro* diagnostic use in calibrating the ADVIA Centaur CoV2Ag assay using the ADVIA Centaur® systems.

Material Description

Material Description	Storage	Stability
ADVIA Centaur CoV2Ag CAL^{a, b}	Frozen at ≤ -20°C	Until expiration date on product
CoV2Ag CAL L: 1.0 mL/vial; frozen Buffer; bovine serum albumin; preservatives	Thawed at 2–8°C	5 days ^c
CoV2Ag CAL H: 1.0 mL/vial; frozen Recombinant SARS-CoV-2 nucleocapsid antigen; buffer; bovine serum albumin; preservatives	At room temperature	4 hours

^a Store in an upright position.

^b Prevent exposure to sunlight and heat.

^c Product may be refrozen 1 time within 5 days for future use. Once re-thawed, use within the same day and discard any remaining material.

Warnings and Precautions

FOR USA:

For Use Under Emergency Use Authorization Only

For *in vitro* diagnostic use.

For Prescription Use.

For Professional Use.

This product has not been FDA cleared or approved but has been authorized for emergency use by FDA under an EUA for use by authorized laboratories; use by laboratories certified under CLIA that meet requirements to perform moderate or high complexity tests.

This product is for use with a test authorized only for the detection of proteins from SARS-CoV-2, not for any other viruses or pathogens.

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of *in vitro* diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb3(b)(1), unless the declaration is terminated, or authorization is revoked sooner.

CAUTION

Federal (USA) law restricts this device to sale by or on the order of a licensed healthcare professional.

Safety data sheets (SDS) available on [siemens-healthineers.com](https://www.siemens-healthineers.com).

CAUTION

This device contains material of animal origin and should be handled as a potential carrier and transmitter of disease.

Dispose of hazardous or biologically contaminated materials according to the practices of your institution. Discard all materials in a safe and acceptable manner and in compliance with prevailing regulatory requirements.

Storage and Stability

Store all calibrators in an upright position, away from light and heat. Do not use products beyond the expiration date printed on the product labeling.

For information about product storage and stability, refer to *Material Description*.

Performing Calibration

Calibration Frequency

For information about calibration frequency, refer to the assay instructions for use.

Preparing the Calibrators

Allow the calibrators to thaw. Gently mix and invert the vials to ensure homogeneity of the material.

Note Use calibrators within the stability limits specified in *Material Description* and discard any remaining material.

Calibration Procedure

The required sample volume for testing depends on several factors. For information about sample volume requirements, refer to the system online help. Perform the calibration procedure using the following steps:

1. Ensure that the calibrator definitions are defined, and that the calibrator values are entered on the system using the assigned value sheet provided.
2. Ensure that the required reagents are loaded for the assay.
3. Schedule the calibrators to the worklist.
4. Label two sample containers with barcode labels: one sample container with a Low Calibrator barcode label, and one sample container with a High Calibrator barcode label.

Note Barcode labels are lot-specific. Do not use barcode labels from one lot of calibrators with any other lot of calibrators.

5. Gently mix the Low and High Calibrators and dispense at least 0.5 mL into the appropriate sample container. Avoid bubbles.

Note This procedure uses sufficient volumes to test each product in triplicate.

6. Load the samples according to the system online help.

Note Dispose of any calibrator remaining in the sample cups after 4 hours. Do not refill calibrator sample cups when the contents are depleted. If required, dispense fresh calibrators.

















Technical Assistance


















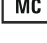


For customer support, contact your local technical support provider or distributor.

siemens-healthineers.com

Definition of Symbols

The following symbols may appear on the product labeling:

Symbol	Symbol Title	Symbol	Symbol Title
	Manufacturer		Authorized representative in the European Community
	Use-by date		Batch code
	Catalog number		Contains sufficient for <n> tests
	Consult Instructions for Use		Version of Instructions for Use
	Internet URL address to access the electronic instructions for use		Revision
	In vitro diagnostic medical device		Unique Device Identifier
	Prescription device (US only)		CE Marking
	CE Marking with Notified Body		Keep away from sunlight


Symbol	Symbol Title	Symbol	Symbol Title
	Temperature limit		Lower limit of temperature
	Upper limit of temperature		Do not freeze
	Do not re-use		This way up
	Recycle		Caution
	Biological risks		Document face up ^a
	Common Units		International System of Units
YYYY-MM-DD	Date format (year-month-day)	YYYY-MM	Date format (year-month)
	Target		Interval
	Handheld barcode scanner		Variable hexadecimal number that ensures the Master Curve and Calibrator definition values entered are valid.
	Lot details		Master Curve definition
	Calibrator lot value		Quality control lot value

^a Indicates Assay-eNote

Legal Information

ADVIA Centaur is a trademark of Siemens Healthineers.

© 2022 Siemens Healthineers. All rights reserved.

 Siemens Healthcare Diagnostics Inc.
511 Benedict Avenue
Tarrytown, NY 10591 USA

Siemens Healthineers Headquarters
Siemens Healthcare GmbH
Henkestraße 127
91052 Erlangen
Germany
Phone: +49 9131 84-0
siemens-healthineers.com

SARS-CoV-2 Antigen Quality Control (CoV2Ag QC)

Current Revision and Date^a	Rev. 03, 2022-05
Product Name	ADVIA Centaur SARS-CoV-2 Antigen Quality Control (CoV2Ag QC)
Abbreviated Product Name	ADVIA Centaur CoV2Ag QC
	2 x 1.6 mL negative quality control, level 1 CONTROL - 1 REF 11207868 2 x 1.6 mL positive quality control, level 2 CONTROL + 2 Quality control assigned value sheet
Systems	ADVIA Centaur systems

^a A vertical bar in the page margin indicates technical content that differs from the previous version.

FOR USA:

For Use Under Emergency Use Authorization Only

For *in vitro* diagnostic use.

For Prescription Use.

For Professional Use.

Intended Use

The ADVIA Centaur® SARS-CoV-2 Antigen Quality Control (CoV2Ag QC) is for *in vitro* diagnostic use in monitoring the performance of the ADVIA Centaur CoV2Ag assay using the ADVIA Centaur® systems for the qualitative detection of SARS-CoV-2 nucleocapsid antigen.

This test has not been FDA cleared or approved; but has been authorized by FDA under an Emergency Use Authorization (EUA) for use by authorized laboratories; use by laboratories certified under CLIA that meet the requirements to perform moderate or high complexity tests.

This test has been authorized only for the detection of proteins from SARS-CoV-2, not for any other viruses or pathogens. The emergency use of this test is authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of *in vitro* diagnostic tests for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Act, 21 U.S.C. § 360bbb3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

Material Description

Material Description	Storage	Stability
ADVIA Centaur CoV2Ag QC^{a, b}	Frozen at ≤ -20°C	Until expiration date on product
CoV2Ag Control 1: 1.6 mL/vial; frozen Buffer; bovine serum albumin; preservatives	Thawed at 2–8°C	7 days ^c
CoV2Ag Control 2: 1.6 mL/vial; frozen Recombinant SARS-CoV-2 nucleocapsid antigen; buffer; bovine serum albumin; preservatives	At room temperature	4 hours

^a Store in an upright position.

^b Prevent exposure to sunlight and heat.

^c Product may be refrozen 1 time within 7 days for future use. Once re-thawed, use within the same day and discard any remaining material.

Warnings and Precautions

FOR USA:

For Use Under Emergency Use Authorization Only

For *in vitro* diagnostic use.

For Prescription Use.

For Professional Use.

This product has not been FDA cleared or approved but has been authorized for emergency use by FDA under an EUA for use by authorized laboratories; use by laboratories certified under CLIA that meet requirements to perform moderate or high complexity tests.

This product is for use with a test authorized only for the detection of proteins from SARS-CoV-2, not for any other viruses or pathogens.

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of *in vitro* diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb3(b)(1), unless the declaration is terminated, or authorization is revoked sooner.

CAUTION

Federal (USA) law restricts this device to sale by or on the order of a licensed healthcare professional.

Safety data sheets (SDS) available on [siemens-healthineers.com](https://www.siemens-healthineers.com).

CAUTION

This device contains material of animal origin and should be handled as a potential carrier and transmitter of disease.

Dispose of hazardous or biologically contaminated materials according to the practices of your institution. Discard all materials in a safe and acceptable manner and in compliance with prevailing regulatory requirements.

Storage and Stability

Store quality control materials in an upright position, away from light and heat. Do not use products beyond the expiration date printed on the product labeling.

For information about product storage and stability, refer to *Material Description*.

Performing Quality Control

Perform the quality control procedure at least once during each day that samples are analyzed. Test quality control samples after a successful calibration.

Follow government regulations or accreditation requirements for quality control frequency. Individual laboratory quality control programs and procedures may require more frequent quality control testing.

Treat all quality control samples the same as patient samples.

Preparing the Quality Control Materials

Allow quality control materials to thaw. Gently mix and invert the vials to ensure homogeneity of the material.

Note Use quality control material within the stability limits specified in *Material Description* and discard any remaining material.

Quality Control Procedure

The required sample volume for testing depends on several factors. For information about sample volume requirements, refer to the system online help. Perform the quality control procedure using the following steps:

1. Ensure that the quality control definitions are defined, and that the quality control values are entered on the system using the assigned value sheet provided.
2. Ensure that the required reagents are loaded for the assay.
3. Schedule the quality control samples to the worklist.
4. Label two sample containers with barcode labels: one sample container for the positive control, and one sample container for the negative control.

Note Barcode labels are lot-specific. Do not use barcode labels from one lot of controls with any other lot of controls.

5. Gently mix each vial of quality control material and dispense at least 0.4 mL into the appropriate sample container. Avoid bubbles.

Note This procedure uses sufficient volumes to test each product in duplicate.

6. Load the samples according to the system online help.

Note Dispose of any QC material that remains in the sample container after 4 hours. Do not refill or reuse sample containers. Do not return any QC material back into the original container.

Taking Corrective Action

If the quality control results do not fall within the assigned values, do not report results. Perform corrective actions in accordance with established laboratory protocol. For suggested protocol, refer to the system online help.

Expected Values

For the assigned values, refer to the quality control assigned value sheet provided. A satisfactory level of performance is achieved when the analyte values obtained are within the expected control interval for the system or within your interval, as determined by an appropriate internal laboratory quality control scheme. Follow your laboratory's quality control procedures if the results obtained do not fall within the acceptable limits. For information about entering QC definitions, refer to the system online help.

The assigned values are traceable to the standardization of the assay. For additional information, refer to the assay instructions for use.

Limitations

The ADVIA Centaur CoV2Ag QC is for use only with the ADVIA Centaur CoV2Ag assay. Assay values have not been established for assays other than the ADVIA Centaur CoV2Ag assay.

The results obtained using quality control material depend on several factors. Erroneous results can occur from causes such as improper storage, inadequate mixing, reconstitution errors, or sample handling errors associated with system or assay procedures.

The assigned control values should be used as a guide in evaluating performance. The control targets and intervals should be adapted to each laboratory's individual requirements. Values obtained should fall within the established interval. Each laboratory should establish corrective measures if individual values fall outside the interval. Follow the applicable government regulations and local guidelines for quality control.























Technical Assistance














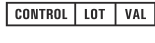
For customer support, contact your local technical support provider or distributor.

siemens-healthineers.com

Definition of Symbols

The following symbols may appear on the product labeling:

Symbol	Symbol Title	Symbol	Symbol Title
	Manufacturer		Authorized representative in the European Community
	Use-by date		Batch code
	Catalog number		Contains sufficient for <n> tests
	Consult Instructions for Use		Version of Instructions for Use
	Internet URL address to access the electronic instructions for use		Revision
	In vitro diagnostic medical device		Unique Device Identifier
	Prescription device (US only)		CE Marking
	CE Marking with Notified Body		Keep away from sunlight
	Temperature limit		Lower limit of temperature
	Upper limit of temperature		Do not freeze
	Do not re-use		This way up

Symbol	Symbol Title	Symbol	Symbol Title
	Recycle		Caution
	Biological risks		Document face up ^a
	Common Units		International System of Units
YYYY-MM-DD	Date format (year-month-day)	YYYY-MM	Date format (year-month)
	Target		Interval
	Handheld barcode scanner		Variable hexadecimal number that ensures the Master Curve and Calibrator definition values entered are valid.
	Lot details		Master Curve definition
	Calibrator lot value		Quality control lot value

^a Indicates Assay-eNote

Legal Information

ADVIA Centaur is a trademark of Siemens Healthineers.

© 2022 Siemens Healthineers. All rights reserved.



Siemens Healthcare Diagnostics Inc.
511 Benedict Avenue
Tarrytown, NY 10591 USA

Siemens Healthineers Headquarters

Siemens Healthcare GmbH
Henkestraße 127
91052 Erlangen
Germany
Phone: +49 9131 84-0
siemens-healthineers.com

IVD

FOR US

RxOnly**For Emergency Use Authorization Only**

- This product has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA for use by authorized laboratories.
- This product has been authorized only for the detection of proteins from SARS-CoV-2, not for any other viruses or pathogens.
- The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of *in vitro* diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3 (b)(1), unless the declaration is terminated or authorization is revoked sooner.
- This card is not the full instructions for use (IFU). The full IFU can be downloaded from the Siemens Healthineers website at [siemens.com/eifu](https://www.siemens.com/eifu); a printed copy of the IFU can be obtained free of charge by contacting Siemens Healthineers Customer Support at 1-888-588-3916.

ADVIA Centaur® SARS-CoV-2 Antigen (CoV2Ag) Calibrator (CAL)
ADVIA Centaur® SARS-CoV-2 Antigen (CoV2Ag) Quality Control (QC)
Sample Inactivation Media

REF 11208047

REF 11207868

REF 11208651

IVD

FOR US

RxOnly

For Emergency Use Authorization Only

- This product has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA for use by authorized laboratories.
- This product is for use with a test authorized only for the detection of proteins from SARS-CoV-2, not for any other viruses or pathogens.
- The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of *in vitro* diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3 (b)(1), unless the declaration is terminated or authorization is revoked sooner.
- This card is not the full instructions for use (IFU). The full IFU can be downloaded from the Siemens Healthineers website at [siemens.com/eifu](https://www.siemens.com/eifu); a printed copy of the IFU can be obtained free of charge by contacting Siemens Healthineers Customer Support at 1-888-588-3916.

Sample Inactivation Media

For Use Under Emergency Use Authorization Only

For *in vitro* diagnostic use.

For prescription use only.

Current Revision and Date ^a	Rev. 02, 2022-05
Product Name	Sample Inactivation Media (50 x 1.0 mL) REF 11208651
Systems	Atellica IM Analyzer ADVIA Centaur XP system ADVIA Centaur XPT system

^a A vertical bar in the page margin indicates technical content that differs from the previous version.

Intended Use

Sample Inactivation Media is used for the processing of anterior nasal swab specimens prior to testing samples with the Atellica® IM SARS-CoV-2 Antigen (CoV2Ag) assay or ADVIA Centaur® SARS-CoV-2 Antigen (CoV2Ag) assay.

Material Description

Material Description	Storage	Stability
Sample Inactivation Media^a 1.0 mL/tube Buffer; surfactant; fetal bovine serum; preservatives	Unopened at 2–30°C	Until expiration date on product

^a Prevent exposure to sunlight and heat.

Warnings and Precautions

For Use Under Emergency Use Authorization Only

For *in vitro* diagnostic use only.

For prescription use only.

This product has not been FDA cleared or approved, but has been authorized by FDA under an EUA for use by authorized laboratories; use by laboratories certified under CLIA that meet requirements to perform moderate or high complexity tests.

This product has been authorized only for the detection of proteins from SARS-CoV-2, not for any other viruses or pathogens.

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of *in vitro* diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

CAUTION

Federal (USA) law restricts this device to sale by or on the order of a licensed healthcare professional.

Safety data sheets (SDS) available on [siemens-healthineers.com](https://www.siemens-healthineers.com).

CAUTION

This device contains material of animal origin and should be handled as a potential carrier and transmitter of disease.

Dispose of hazardous or biologically contaminated materials according to the practices of your institution. Discard all materials in a safe and acceptable manner and in compliance with prevailing regulatory requirements.

Storage and Stability

Store Sample Inactivation Media at 2–30°C, away from light and heat.

Avoid bacterial contamination.

Do not use product beyond the expiration date printed on the product labeling.

For information about product storage and stability, refer to *Material Description*.

Sample Inactivation Procedure

For information about using the Sample Inactivation Media to process samples prior to testing, refer to the assay instructions for use.

Limitations

- The performance of the Sample Inactivation Media has not been established with any assays or specimen types other than those defined in the *Intended Use*.
- Do not use to inactivate any viruses other than SARS-CoV-2.
- All samples (including inactivated samples) and reagents containing biological materials used for the assay must be considered as potentially able to transmit infectious agents. All samples, reagents, and associated waste must be handled with extreme care and disposed of in compliance with the laboratory guidelines and the statutory provisions in force in each country.







































Technical Assistance

For customer support, contact your local technical support provider or distributor.

[siemens-healthineers.com](https://www.siemens-healthineers.com)

Definition of Symbols

The following symbols may appear on the product labeling:

Symbol	Symbol Title	Symbol	Symbol Title
	Manufacturer		Authorized representative in the European Community
	Use-by date		Batch code
	Catalog number		Contains sufficient for <n> tests
	Consult Instructions for Use		Version of Instructions for Use
	Internet URL address to access the electronic instructions for use		Revision
	<i>In vitro</i> diagnostic medical device		Unique Device Identifier
	Prescription device (US only)		CE Marking
	CE Marking with Notified Body		Keep away from sunlight
	Temperature limit		Lower limit of temperature
	Upper limit of temperature		Do not freeze
	Do not re-use		This way up
	Recycle		Caution
	Biological risks		Document face up ^a
	Common Units		International System of Units
	Date format (year-month-day)		Date format (year-month)
	Handheld barcode scanner		Mixing of substances
	Target		Interval
	Variable hexadecimal number that ensures the Master Curve and Calibrator definition values entered are valid.		Material
	Unique material identification number		Name of control

Symbol	Symbol Title	Symbol	Symbol Title
CONTROL TYPE	Type of control	CAL LOT VAL	Calibrator lot value
CONTROL LOT VAL	Quality control lot value		

^a Indicates Assay-eNote

Legal Information

Atellica and ADVIA Centaur are trademarks of Siemens Healthineers.

© 2022 Siemens Healthineers. All rights reserved.



Siemens Healthcare Diagnostics Inc.
511 Benedict Avenue
Tarrytown, NY 10591 USA

Siemens Healthineers Headquarters

Siemens Healthcare GmbH
Henkestraße 127
91052 Erlangen
Germany
Phone: +49 9131 84-0
siemens-healthineers.com